

LOOK™ SUTURE

# PRODUCT CATALOG



- Consistency
- Quality
- Value

# WELCOME

---

The LOOK™ line of sutures is one of the best known and most trusted suture brands on the market. From the beginning, these sutures have been known for their quality—and that legacy of quality continues to meet the demands of surgery professionals worldwide. The precision manufacturing of our LOOK™ brand sutures ensures consistent needle strength and durability that holds up to every pass.

LOOK™ Sutures are offered in a variety of absorbable and non-absorbable materials and are attached to needles with superior sharpness, penetration, and strength characteristics.

Surgical Specialties Corporation is proud to provide the LOOK™ suture, along with flexible and accommodating customer service, and value to its customers.

# Table Of Contents

---

How to Use Our Catalog .....	1
Materials Chart .....	2
A Guide To Surgical Sutures .....	3
Needle Anatomy and Terminology .....	4
Needle Shapes And General Uses .....	5
Needle Points and Cutting Edges .....	6
Needle Charts .....	7
<b>Absorbable Sutures</b>	
Plain Gut.....	9
Chromic Gut .....	11
PolySyn FA™ (polyglycolic acid) .....	13
PolySyn™ (polyglycolic acid) .....	14
<b>Non-Absorbable Sutures</b>	
Silk .....	18
Nylon.....	21
Polypropylene.....	24
Polyviolene™ (Polyester) - Uncoated.....	26
Additional Products .....	27
Product Index.....	31

# Introduction

## How to Use Our Catalog

The pages in this catalog show the following information for each product code. Some products feature an expanded chart with detailed measurements (not shown in this sample) to ensure proper needle selection.

**Needle Description** - Contains an illustration of the needle. The LOOK™ needle code is purple, the comparable Ethicon™ needle code is black, and the comparable Syneture™ needle code is dark green.


**Needle Geometry** - Type of needle with tip geometry

**Material** - Color and style of suture

Precision  
Reverse Cutting



**Length** - Length in inches/centimeters

Needle	Material	Length	Size				
			5-0	4-0	3-0	2-0	0
<b>PC31</b> 19mm 3/8 Circle  PS-2 P-12	Undyed Braid	18"/45cm	495B J495G SL5636G	496B • SL5627G			

**Size Grid** - Indicates the U.S.P. size for each product code listed.

• • indicates that other manufacturers do not have a comparable product code or needle. LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

LOOK™ product code

Ethicon™ comparable product code




Syneture™ comparable product code

# Materials Chart






LOOK™ Suture	Ethicon™	Syneture™	Material
Plain Gut	Surgical Gut	Plain Gut	Natural Absorbable Gut
Chromic Gut	Surgical Gut	Chromic Gut	Natural Absorbable Gut
PolySyn FA™	Vicryl Rapide™	Velosorb™ Fast	Synthetic Absorbable, Braided, Coated Polyglycolic Acid
PolySyn™	Vicryl™	Polysorb™/Dexon™	Synthetic Absorbable, Braided, Coated Polyglycolic Acid
Silk	Permahand Silk™	SofSilk™	Non-Absorbable, Braided, Coated Silk
Nylon	Ethilon™	Monosof™/Dermalon™	Non-Absorbable, Monofilament Nylon
Polypropylene	Prolene™	Surgipro™ II	Non-Absorbable, Monofilament Polypropylene
Uncoated Polyviolene™	Mersilene™	Surgidac™	Non-Absorbable, Braided, Uncoated Polyester
Coated Polyviolene™	Ethibond™	Ti-Cron™	Non-Absorbable, Braided, Coated Polyester

# A Guide To Surgical Sutures

## Absorbable


Material	Material Description	BSR Profile	Mass Absorption
 <p><b>Gut</b></p> <ul style="list-style-type: none"> <li>• Surgical Gut Suture - Plain/Chromic</li> <li>* Plain Gut/Chromic Gut</li> </ul>	<p><b>Twisted</b></p> <ul style="list-style-type: none"> <li>• Plain &amp; Chromic</li> <li>• Packaged Dry or in Wetting Solution</li> </ul> <p><b>Common uses:</b> General Closure, OB/Gyn, Urology, Bowel Anastomosis</p>	<p><b>In vivo strength retention:</b> Varies</p>	<p>Varies depending upon type of suture (plain vs. chromic) and site of placement</p>
 <p><b>PolySyn FA™</b> (glycolic acid)</p> <ul style="list-style-type: none"> <li>• VICRYL RAPIDE™ * Velosorb™ Fast</li> </ul>	<p><b>Braided</b></p> <ul style="list-style-type: none"> <li>• Undyed (clear) • Coated</li> </ul> <p><b>Common uses:</b> General Closure, OB/Gyn, Urology</p>	<p><b>In vivo strength retention:</b> 48% at 7 days 0% at 14 days</p>	<p>Essentially complete between 49 and 63 days</p>
 <p><b>PolySyn™</b> (glycolic acid)</p> <ul style="list-style-type: none"> <li>• VICRYL™ * Polysorb™ / Dexon™</li> </ul>	<p><b>Braided</b></p> <ul style="list-style-type: none"> <li>• Undyed (clear) • Dyed (violet)</li> </ul> <p><b>Monofilament</b></p> <ul style="list-style-type: none"> <li>• Undyed (clear) • Dyed (violet)</li> <li>• Coated &amp; Uncoated</li> </ul> <p><b>Common uses:</b> Skin Closure, OB/Gyn, Urology, Gastrointestinal</p>	<p><b>In vivo strength retention:</b> 50% at 14 days 20% at 21 days</p>	<p>Essentially complete between 50 and 90 days</p>

## Non-Absorbable

Material	Material Description	Qualities
 <p><b>Silk</b></p> <ul style="list-style-type: none"> <li>• PERMA-HAND™ Silk * SofSilk™</li> </ul>	<p><b>Braided</b></p> <ul style="list-style-type: none"> <li>• Undyed (white) • Dyed (black) • Coated</li> </ul> <p><b>Common uses:</b> General &amp; Skin Closure, Gastrointestinal, Cardiovascular Surgery, Plastic Surgery, Ophthalmology</p>	<ul style="list-style-type: none"> <li>• Excellent handling and tying characteristics</li> </ul>
 <p><b>Nylon</b></p> <ul style="list-style-type: none"> <li>• ETHILON™ * Monosof™/Dermalon™</li> </ul>	<p><b>Monofilament</b></p> <ul style="list-style-type: none"> <li>• Undyed (clear) • Dyed (blue) (black)</li> <li>• Packaged dry or in wetting solution</li> </ul> <p><b>Common uses:</b> General &amp; Skin Closure, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery, Gastrointestinal</p>	<ul style="list-style-type: none"> <li>• Passes easily through tissue</li> </ul>
 <p><b>Polypropylene</b></p> <ul style="list-style-type: none"> <li>• PROLENE™ * Surgipro™/Surgipro™ II</li> </ul>	<p><b>Monofilament</b></p> <ul style="list-style-type: none"> <li>• Undyed (clear) • Dyed (blue)</li> </ul> <p><b>Common uses:</b> General &amp; Skin Closure, Neurosurgery, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery</p>	<ul style="list-style-type: none"> <li>• Inert</li> <li>• Easy to handle</li> <li>• Passes easily through tissue</li> </ul>
 <p><b>Polyviolene™ (Uncoated)</b> (poly (ethylene-terephthate))</p> <ul style="list-style-type: none"> <li>• MERSILENE™ * Surgidac™</li> </ul>	<p><b>Braided</b></p> <ul style="list-style-type: none"> <li>• Undyed (white) • Dyed (green)</li> <li>• Uncoated</li> </ul> <p><b>Common uses:</b> General &amp; Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures</p>	<ul style="list-style-type: none"> <li>• Smooth</li> <li>• Easy to handle</li> <li>• Good knot security</li> </ul>
 <p><b>Polyviolene™ (Coated)</b> (poly (ethylene-terephthate))</p> <ul style="list-style-type: none"> <li>• ETHIBOND EXCEL™ * Ti-Cron™</li> </ul>	<p><b>Braided</b></p> <ul style="list-style-type: none"> <li>• Undyed (white) • Dyed (green)</li> <li>• Coated</li> </ul> <p><b>Common uses:</b> General &amp; Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures</p>	<ul style="list-style-type: none"> <li>• Smooth</li> <li>• Easy to handle</li> <li>• Good knot security</li> </ul>

• Comparable Ethicon™ Material \* Comparable Synture™ Material

# Needle Anatomy And Terminology

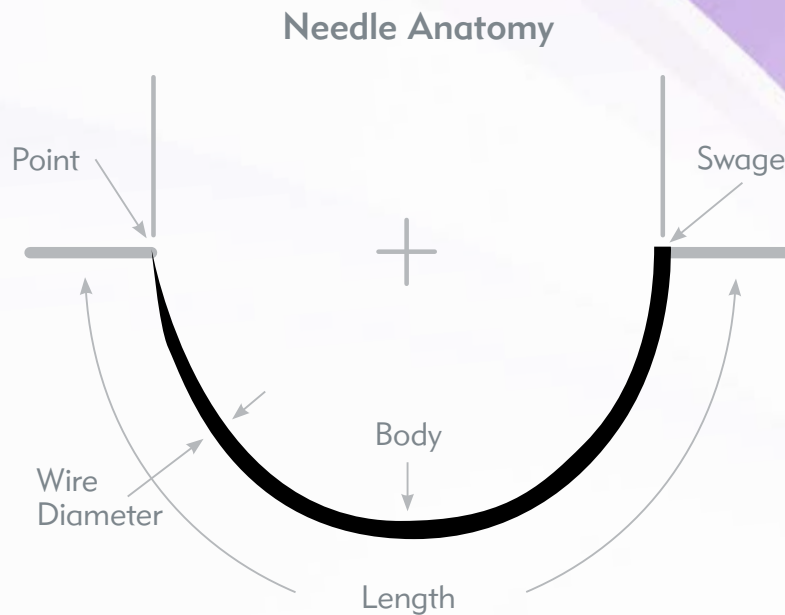

## Needle Terminology

**Length:** Distance measured along the circumference of the needle, from the point to the swage. It is depicted on the packaging.

**Shape:** Different surgeries or surgical techniques demand different needle shapes.

**Point:** The sharpened part of the needle that first penetrates the tissue. The point consists of the area from the sharpened tip of the needle to the widest cross section of the needle body. The point determines ease of penetration and the initial size and shape of the hole.

**Swage:** Connection point of the suture and the needle. Micro teeth grip the suture material and prevent it from slipping out.

**Wire Diameter:** Gauge or thickness of the wire originally used to form the needle body.

**Body:** The part between the point and the swage of the needle.


**Included Angle (Curvature):** The fraction of the circle included in the curve of the needle. This is depicted on the packaging.

# Needle Shapes And General Uses


Straight (AKA: KEITH): Skin Closure

3/8 Circle: Skin Closure, Ophthalmic


1/2 Circle: Skin Closure, GI, Cardiovascular, Ophthalmic

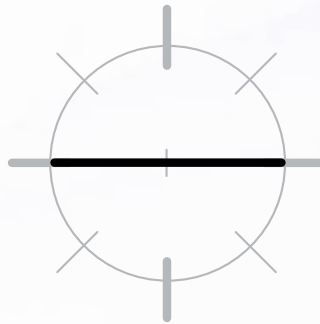

Straight

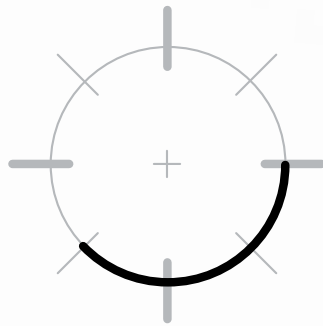



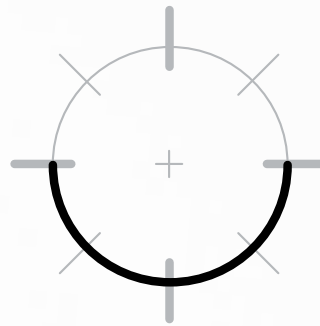





3/8 Circle



1/2 Circle



# Needle Points And Cutting Edges

Needle Type	Point Geometry	Description	Commonly Used For
Conventional Cutting		Triangular point with cutting edge on the inner curvature	General skin closure, subcutaneous tissue, sometimes for ophthalmic surgery, plastic or reconstructive surgery
Precision Conventional Cutting		Triangular point with cutting edge on the inner curvature and precision-honed for added sharpness	General skin closure, subcutaneous tissue, sometimes for ophthalmic surgery, plastic or reconstructive surgery
Reverse Cutting		Triangular point with cutting edge on the outer curvature	Skin closure, ligamentous or fibrous tissues
Precision Reverse Cutting		Triangular point with cutting edge on the outer curvature and precision-honed for added sharpness	Primarily used in plastic surgery for delicate procedures and where a good cosmetic result is desired
Taper Point		Round shaft, straight or curved, taper point with no cutting edge	Soft tissue closure such as fascia, vascular, gastrointestinal, and most soft tissue below the skin
Taper Cutting		Four (4) cutting edges with a taper shaft	For penetration through tough tissue
Side Cutting Lancet		Spatula needle with the cutting edge on the inner, or outer curvature	Ophthalmic and microsurgery
Straight Cutting		Reverse cutting tip on a straight shaft	General skin closure

# Needle Chart

## Reverse Cutting Needles



LOOK™	Ethicon™	Syneture™
<b>3/8 Circle</b>		
C16	C-2	C-1 ————— 11mm
C17	C-2	C-1 ————— 12mm
•	C-3	C-12 ————— 14mm
C22	FS-3	C-12 ————— 16mm
C6	FS-2	C-13 ————— 18mm
C7	FS-1	C-14 ————— 24mm
C8	FS	C-15 ————— 26mm
C10	FSL	C-16 ————— 30mm
<b>1/2 Circle</b>		
C26	M-1	C-21 ————— 15mm
C19	•	• ————— 18mm
C20	J-1	C-22 ————— 19mm
C31	X-1	C-23 ————— 24mm
C12	CP-1	GS-11 ————— 27mm
C21	CP-2	GS-10 ————— 26mm
C13	CP-1	GS-11 ————— 26mm
C14	CP	GS-12 ————— 37mm
		————— 40mm

## Precision Reverse Cutting Needles



LOOK™	Ethicon™	Syneture™
<b>3/8 Circle</b>		
C1	P-1	P-10 ————— 11mm
C3	P-3	P-13 ————— 13mm
PC34	PS-3	P-11 ————— 16mm
PC31	PS-2	P-12 ————— 18mm
PC33	PS-1	P-14 ————— 24mm

## Precision Conventional Cutting Needles



LOOK™	Ethicon™	Syneture™
<b>3/8 Circle</b>		
CP1	PC-1	P-13 ————— 11mm

# Needle Chart

## Taper Point Needles



LOOK™	Ethicon™	Syneture™
1/2 Circle		
T22	TF	CV-22
T28	RB-1	CV-23
T14	SH-1	CV-24
T4	SH-1	CV-25
T	SH	V-20
T5	CT-2	GS-22
HR30	V-30	•
T9	MH	V-26
T8	CT-1	GS-21
T10	CT	GS-24

## Straight Cutting Needles



LOOK™	Ethicon™	Syneture™
SC2	KS	SC-2

## Straight Taper Point



LOOK™	Ethicon™	Syneture™
S	•	•

## Taper Cutting Needles



LOOK™	Ethicon™	Syneture™
1/2 Circle		
Y5	V-5	KV-5
HRT31	V-34	KV-34
HRT36	V-34	KV-34

## Side Cutting Lancet





LOOK™	Ethicon™	Syneture™
VSP8	S-24	SS-24

# Plain Gut


## Precision Conventional Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	1
<b>CP1</b>  11mm 3/8 Circle PC-1 P-13	Plain Gut	18"/45cm	<b>1415B</b> 						


## Precision Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C3</b>  13mm 3/8 Circle P-3 P-13	Plain Gut	18"/45cm		<b>596B</b> 686G SG5686	<b>1241B</b> 1644G SG5654G			

## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C17</b>  12mm 3/8 Circle C-2 C-1	Plain Gut	10"/25cm		<b>521B</b> 	<b>531B</b> 			
		18"/45cm		<b>535B</b> 686G CG172	<b>536B</b> 1644G SG721			
<b>C6</b>  18mm 3/8 Circle FS-2 C-13	Plain Gut	18"/45cm		<b>500B</b>  SG751	<b>591B</b>  SG752	<b>592B</b>  SG5630G		
		27"/70cm		<b>550B</b> H820G SG820G	<b>551B</b> H821H SG821G	<b>554B</b> H822H SG822G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).





• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Plain Gut

Reverse Cutting

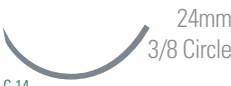

Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C7</b>  24mm 3/8 Circle FS-1 C-14	Plain Gut	27"/70cm				<b>553B</b> H810H SG810	<b>555B</b> H823H SG823	
<b>C26</b>  15mm 1/2 Circle M-1 C-21	Plain Gut	18"/45cm			<b>537B</b> • •			
<b>C31</b>  24mm 1/2 Circle X-1 C-23	Plain Gut	18"/45cm			<b>599B</b> • SG1633G	<b>594B</b> 612G SG612G		
<b>C21</b>  26mm 1/2 Circle CP-2 GS-10	Plain Gut	18"/45cm				<b>593B</b> 612G SG612G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Chromic Gut

## Precision Reverse Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
<b>C1</b>  11mm 3/8 Circle P-1 P-10	Chromic Gut	12"/30cm	<b>1242B</b> 796G G1780						
<b>C3</b>  13mm 3/8 Circle P-3 P-13	Chromic Gut	18"/45cm		<b>1248B</b> 687G SG5687G	<b>595B</b> 1654G SG5644G				
<b>PC31</b>  18mm 3/8 Circle PS-2 P-12	Chromic Gut	18"/45cm			<b>1253B</b> 1637G SG5637G				

## Reverse Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
<b>C16</b>  11mm 3/8 Circle C-2 C-1	Chromic Gut	18"/45cm		<b>543B</b> 792G G1792K					
<b>C17</b>  12mm 3/8 Circle C-2 C-1	Chromic Gut	18"/45cm		<b>563B</b> 687G SG197	<b>546B</b> 1654G SG5644G				
<b>C6</b>  18mm 3/8 Circle FS-2 C-13	Chromic Gut	18"/45cm			<b>558B</b> 1637G SG675	<b>525B</b> 1638H SG676			
		27"/70cm		<b>557B</b> 634G SG634G	<b>559B</b> 635H SG635G	<b>560B</b> 636H SG636G			
<b>C7</b>  24mm 3/8 Circle FS-1 C-14	Chromic Gut	27"/70cm					<b>561B</b> 637H SG637G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).




• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Chromic Gut



## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C8</b>  26mm 3/8 Circle FS C-15	Chromic Gut	27"/70cm						562B • •
<b>C26</b>  15mm 1/2 Circle M-1 C-21	Chromic Gut	18"/45cm		598B • •	544B 1653G •			
<b>C20</b>  20mm 1/2 Circle J-1 C-22	Chromic Gut	18"/45cm			547B 1643G •			

## Taper Point



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>T28</b>  18mm 1/2 Circle RB-1 CV-23	Chromic Gut	14"/35cm		590B • UG202				
<b>T4</b>  22mm 1/2 Circle SH-1 CV-25	Chromic Gut	14"/35cm			589B G181H GG181			

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# PolySyn FA™ / Polyglycolic Acid


## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C6</b>  18mm 3/8 Circle FS-2 C-13	Undyed Braided	27"/70cm				MC-423NC • •		
<b>C7</b>  24mm 3/8 Circle FS-1 C-14	Undyed Braided	27"/70cm				MC-442NC VC442H SL693		

## Taper Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>HRT 36</b>  36mm 1/2 Circle V-34 KV-34	Undyed Braided	36"/90cm					MC-917NC VR917 J9962T	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.



# PolySyn™/ Polyglycolic Acid

## Precision Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C3</b>  13mm 3/8 Circle P-3 P-13	Undyed Braided	18"/45cm	<b>492B</b> J492G SL5688G	<b>493B</b> J493G SL5687G	<b>494B</b> J494G SL5691G			
<b>PC31</b>  19mm 3/8 Circle PS-2 P-12	Undyed Braided	18"/45cm		<b>495B</b> J495G SL5636G	<b>496B</b> J496G SL5627G	<b>497B</b> J497G SL5628G		
<b>PC33</b>  24mm 3/8 Circle PS-1 P-14	Undyed Braided	18"/45cm				<b>683B</b> J683G SL5640G		

## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C17</b>  12mm 3/8 Circle C-2 C-1	Undyed Braided	10"/25cm		<b>385B</b> VC493G SL606	<b>386B</b> VC494G SL605			
<b>C6</b>  19mm 3/8 Circle FS-2 C-13	Violet Braided	27"/70cm			<b>MC-397NC</b> • •			
	Undyed Braided	27"/70cm		<b>421B</b> J421H SL690G	<b>422B</b> <b>422BC</b> J422H SL691G	<b>423B</b> J423H SL822G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# PolySyn™ / Polyglycolic Acid


## Reverse Cutting



Needle	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
<b>C7</b> 24mm 3/8 Circle  FS-1 C-14	Violet Braided	27"/70cm		GC-451NC • •				
	Undyed Braided	27"/70cm		441B J441H SL692G	GC-452NC 442B J442H SL693	443B J443H SL823		
<b>C14</b> 40mm 1/2 Circle  CP GS-12	Violet Braided	27"/70cm						480B J480H CL805

## Reverse Cutting / Taper Point



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C12/T8</b> 27mm /36.6mm 1/2 Circle  CP-1 GS-11 CT-1 GS-21	Violet Braided	54"/135cm						GC-532NC J532H •

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# PolySyn™/ Polyglycolic Acid

## Taper Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
<b>HRT 31</b>  31mm 1/2 Circle V-34 KV-34	Violet Braided	30"/75cm							<b>GC-936NC</b> W9391 

## Taper Point



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
<b>T22</b>  13mm 1/2 Circle TF CV-22	Undyed Braided	27"/70cm		<b>GC-433NC</b>  					
	Violet Braided	18"/45cm	<b>GC-981NC</b>  						
<b>T28</b>  18mm 1/2 Circle RB-1 CV-23	Violet Braided	30"/75cm			<b>GC-304NC</b>  				







Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# PolySyn™/ Polyglycolic Acid

Taper Point 

Needle	Material	Length	Length					
			5-0	4-0	3-0	2-0	0	1
<b>T14</b>  20mm 1/2 Circle SH-1 CV-24	Violet Braided	27"/70cm		<b>310B/ 310BC</b> J310H GL181	<b>311B/ 311BC</b> J311H GL182			
<b>T5</b>  26.50mm 1/2 Circle CT-2 GS-22	Violet Braided	27"/70cm			<b>332B</b> J332H CL882	<b>333B</b> J333H CL883		
<b>T</b>  26mm 1/2 Circle SH V-20	Violet Braided	27"/70cm			<b>316B/ 316BC</b> J316H GL122	<b>317B/ 317BC</b> J317H GL123		
<b>HR30</b>  30mm 1/2 Circle V-30	Violet Braided	30"/75cm			<b>GC-322BC</b> • •	<b>GC-323NC</b> • •		
<b>T8</b>  36.60mm 1/2 Circle	Violet Braided	27"/70cm				<b>339B</b> J339H CL930	<b>340B</b> J340H CL812	<b>341B</b> J341H CL813
		36"/90cm				<b>GC-354NC</b> • •	<b>GC-518NC</b> • •	
<b>T10</b>  40mm 1/2 Circle CT GS-24	Violet Braided	27"/70cm					<b>352B</b> J352H CL802	<b>353B</b> J353H CL803
		36"/90cm					<b>GC-358NC</b> • •	<b>GC-359NC</b> • •

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).




• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Silk



## Precision Reverse Cutting



Needle	Material	Length	Size						
			8-0	7-0	6-0	5-0	4-0	3-0	
<b>C1</b>  11mm 3/8 Circle P-1 P-10	Black Braided	18"/45cm			<b>1260B</b> 639G SS1639G				
<b>C3</b>  13mm 3/8 Circle P-3 P-13	Black Braided	18"/45cm				<b>1265B</b> 640G SS5640G	<b>792B</b> 641G SS5641G		
<b>PC31</b>  18mm 3/8 Circle PS-2 P-12	Black Braided	18"/45cm					<b>1270B</b> 1677G SS5677G		

## Reverse Cutting



Needle	Material	Length	Size						
			8-0	7-0	6-0	5-0	4-0	3-0	
<b>C16</b>  11mm 3/8 Circle C-2 C-1	Black Braided	18"/45cm							<b>795B</b> • •
<b>C17</b>  12mm 3/8 Circle C-2 C-1	Black Braided	18"/45cm			<b>752B</b> • S1172	<b>774B</b> • S1173	<b>787B</b> 739G S1174		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Silk

## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C22</b> 16mm 3/8 Circle  FS-3 C-12	Black Braided	18"/45cm	<b>750B</b> • SS681					
<b>C6</b> 18mm 3/8 Circle  FS-2 C-13	Black Braided	18"/45cm		<b>754B</b> 682G SS682G	<b>781B</b> 683G SS683G	<b>784B</b> 1679H SS694		
		30"/75cm				<b>768B</b> 1679H SS745		
<b>C7</b> 24mm 3/8 Circle  FS-1 C-14	Black Braided	18"/45cm			<b>776B</b> 629G SS629	<b>786B</b> 684G SS684G		
		30"/75cm				<b>775B</b> • •		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).








• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Silk
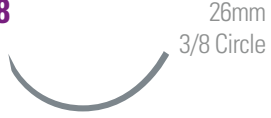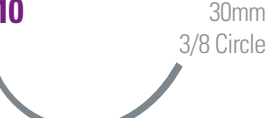
## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C8</b>  26mm 3/8 Circle FS C-15	Black Braided	18"/45cm					<b>770B</b> 685G SS685G	
<b>C10</b>  30mm 3/8 Circle FSL C-16	Black Braided	18"/45cm				<b>765B</b> 675H SS675	<b>772B</b> • •	
<b>C26</b>  15mm 1/2 Circle M-1 C-21	Black Braided	18"/45cm		<b>779B</b> • •	<b>777B</b> M154T •	<b>778B</b> • •		
<b>C19</b>  18mm 1/2 Circle	Black Braided	18"/45cm				<b>783B</b> • •		
<b>C20</b>  20mm 1/2 Circle J-1 C-22	Black Braided	18"/45cm			<b>780B</b> 734G SS733			
<b>C21</b>  26mm 1/2 Circle CP-2 GS-10	Black Braided	18"/45cm				<b>782B</b> • •		
<b>C31</b>  24mm 1/2 Circle X-1 C-23	Black Braided	18"/45cm				<b>785B</b> 632G SS732		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).


• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Nylon






## Precision Conventional Cutting



Needle	Material	Length	Size						
			6-0	5-0	4-0	3-0	2-0	0	
<b>CP1</b>  11mm 3/8 Circle PC-1 P-13	Black Monofilament	18"/45cm	<b>1456B</b> 1856G SN1956		<b>1454B</b> 1854G SN1954				

## Precision Reverse Cutting



Needle	Material	Length	Size						
			7-0	6-0	5-0	4-0	3-0	2-0	
<b>C1</b>  11mm 3/8 Circle P-1 P-10	Black Monofilament	18"/45cm	<b>1276B</b> 1696G SN3686	<b>1279B</b> 697G SN3697					
<b>C3</b>  13mm 3/8 Circle P-3 P-13	Black Monofilament	8"/20cm		<b>1273B</b> • •	<b>1284B</b> • •				
		18"/45cm		<b>1282B</b> 1698G SN5696G	<b>1285B</b> 698G SN5698G	<b>1291B</b> 699G SN5699G			
	Clear Monofilament	18"/45cm			<b>1275B</b> 690G SN5690G				
<b>PC34</b>  16mm 3/8 Circle PS-3 P-11	Black Monofilament	18"/45cm		<b>1280B</b> 1665G SN5767G	<b>1283B</b> 1668G SN5668G				
<b>PC31</b>  18mm 3/8 Circle PS-1 P-14	Black Monofilament	18"/45cm		<b>1281B</b> • SN5665G	<b>1288B</b> 1666G SN5666G	<b>1290B</b> 1667G SN5667G			
<b>PC33</b>  24mm 3/8 Circle PS-1 P-14	Black Monofilament	18"/45cm						<b>1294B</b> 1669G SN5669G	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.







LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.



# Nylon

## Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C17</b> C-2 C-1  12mm 3/8 Circle	Black Monofilament	10"/25cm	<b>911B</b> • SN667G	<b>913B</b> • SN871	<b>931B</b> • SN772G			
	Clear Monofilament	10"/25cm		<b>946B</b> 668G SN5690G				
	Black Monofilament	18"/45cm	<b>916B</b> 667G SN667G	<b>917B</b> 668G SN871	<b>921B</b> 1954G SN872			
<b>C22</b> FS-3 C-12  16mm 3/8 Circle	Black Monofilament	18"/45cm	<b>915B</b> 1665G SN660G	<b>919B</b> 1668G SN659G				
<b>C6</b> FS-2 C-13  18mm 3/8 Circle	Black Monofilament	10"/25cm		<b>945B</b> 661H SN661G	<b>944B</b> 662G SN662G			
		18"/45cm		<b>918B</b> 661G SN661G	<b>922B</b> 662G SN662G	<b>928B</b> • SN653		
<b>C7</b> FS-1 C-14  24mm 3/8 Circle	Black Monofilament	10"/25cm			<b>947B</b> 1629H SN658			
		18"/45cm			<b>920B</b> 1629H SN658	<b>925B</b> 663G SN663G		
<b>C8</b> FS C-15  26mm 3/8 Circle	Black Monofilament	18"/45cm					<b>927B</b> 664G SN764G	
<b>C10</b> FSL C-16  30mm 3/8 Circle	Clear Monofilament	27"/70cm				<b>950B</b> 1671H SN671		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).


• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Nylon


## Straight Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>SC2</b> 60mm Straight  KS SC-2	Black Monofilament	30"/75cm				<b>923B</b> 627H SN627		

## Straight Taper Point



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>S</b> 61mm Straight 	Black Monofilament	30"/75cm					<b>933B</b> • •	


Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.



LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Polypropylene

## Precision Conventional Cutting

Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>CP1</b>  11mm 3/8 Circle PC-1 P-13	Blue Monofilament	8"/20cm	<b>1057B</b> • •	<b>1058B</b> • •				

## Precision Reverse Cutting

Needle	Material	Length	Size					
			7-0	6-0	5-0	4-0	3-0	2-0
<b>C1</b>  11mm 3/8 Circle P-1 P-10	Blue Monofilament	18"/45cm		<b>1079B</b> 8697G SP5697G				
<b>C3</b>  13mm 3/8 Circle P-3 P-13	Blue Monofilament	8"/20cm		<b>1077B</b> • •	<b>1078B</b> • •			
		18"/45cm		<b>1082B</b> 8695G SP5695G	<b>1085B</b> 8698G SP5698G	<b>1091B</b> 8699G SP5699G		
<b>PC31</b>  18mm 3/8 Circle PS-2 P-12	Blue Monofilament	18"/45cm				<b>1020B</b> 8682G SP5682G		

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).


• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Polypropylene

Reverse Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>C17</b>  12mm 3/8 Circle C-2 C-1	Blue Monofilament	10"/25cm	1011B ● ●	1013B ● ●	1031B ● ●			
<b>C6</b>  18mm 3/8 Circle FS-2 C-13	Blue Monofilament	18"/45cm		1018B 8661G SP661G	1022B 8683G SP683G			
<b>C7</b>  24mm 3/8 Circle FS-1 C-14	Blue Monofilament	18"/45cm			1095B ● ●			

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

● indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Polyviolene™ / Polyester - Uncoated


## Reverse Cutting



Needle	Material	Length	Size					
			9-0	8-0	7-0	6-0	5-0	4-0
<b>C22</b>  16mm 3/8 Circle FS-3 C-12	Green Braided Uncoated	18"/45cm					<b>1176B</b> X681G •	
<b>C6</b>  18mm 3/8 Circle FS-2 C-13	Green Braided Uncoated	18"/45cm						<b>1180B</b> X683G •
	White Braided Uncoated	18"/45cm						<b>1205B</b> R633G XTC692


## Side Cutting Lancet



Needle	Material	Length	Size					
			9-0	8-0	7-0	6-0	5-0	4-0
<b>VSP8</b>  8mm 1/4 Circle S-24 SS-24	White Braided Uncoated	18"/45cm					<b>FA-893NC</b> W893 •	

## Taper Cutting



Needle	Material	Length	Size					
			6-0	5-0	4-0	3-0	2-0	0
<b>Y5</b>  16mm 1/2 Circle V-5 KV-5	Green Braided Uncoated	18"/45cm		<b>797B</b> X951H •				

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Additional Products


## Silk Bulk Suture - Reverse Cutting



Needle	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
<b>C7</b>  24mm 3/8 Circle FS-1 C-14	Silk, Black Braided Strands	250 Pieces 18"/45cm			OF766 • •			
<b>C8</b>  26mm 3/8 Circle FS C-15	Silk, Black Braided Strands	250 Pieces 18"/45cm				OF770 • •		
<b>C10</b>  30mm 3/8 Circle FSL C-16	Silk, Black Braided Strands	250 Pieces 18"/45cm					OF772 • •	
<b>C21</b>  26mm 1/2 Circle CP-2 GS-10	Silk, Black Braided Strands	250 Pieces 18"/45cm			OF782 • •			

## Silk Bulk Suture - Straight Cutting



NEEDLE	MATERIAL	LENGTH	SIZE					
			6-0	5-0	4-0	3-0	2-0	0
<b>SC2</b>  60mm Straight KS SC-2	Silk, Black Braided Strands	250 Pieces 30"/75cm				OF762 • •	OF767 • •	

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).



• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Additional Products

## Nylon Bulk Suture - Reverse Cutting



Needle	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
<b>C6</b>  FS-2 C-13	Nylon, Black Monofilament Strands	250 Pieces 18"/45cm			<b>OF928</b> • •			
<b>C8</b>  FS C-15	Nylon, Black Monofilament Strands	250 Pieces 18"/45cm				<b>OF927</b> • •		

## LOOK™ Nylon Spools / Non-Sterile

Material	Length	Size					
		5-0	4-0	3-0	2-0	0	1
Non-Absorbable Nylon, Black Braided	100 yards	<b>SP175</b> • •	<b>SP176</b> • •				
Non-Absorbable Nylon, Black Monofilament	100 yards	<b>SP185</b> • •	<b>SP186</b> • •	<b>SP187</b> • •	<b>SP188</b> • •	<b>SP189</b> • •	<b>SP190</b> • •

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Additional Products

## LOOK™ Spools Polyviolene™/ Non-Sterile

Material	Length	Size					
		7-0	6-0	5-0	4-0	3-0	2-0
Non-Absorbable Polypropylene, Blue Monofilament	100 yards			SP135 • •	SP136 • •	SP137 • •	SP138 • •

## LOOK™ Silk Spools / Non-Sterile

Material	Length	Size							
		6-0	5-0	4-0	3-0	2-0	0	1	2
Silk, Black Braided on Spools	25 yards	SP102 • •	SP103 • •	SP104 • •	SP105 • •	SP106 • •			
	50 yards								SP121 • •
	100 yards	SP114 • •	SP115 • •	SP116 • •	SP117 • •	SP118 • •	SP119 • •	SP120 • •	
Silk, White Braided on Spools	2000 yards				SP1230 • •				

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

• • indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.



# Additional Products

## LOOK™ Polyviolene™ Spools / Non-Sterile

Material	Length	Size					
		5-0	4-0	3-0	2-0	0	1
Non-Absorbable Polyester, Green Braided, Coated	250 yards		SP252/ SP252R ● ●				
	3281 yards		SP2520 ● ●				

## LOOK™ Polyviolene™ Spools / Non-Sterile

Material	Length	Size					
		0	1	2	3	4	5
Non-Absorbable Polyester, Green Braided, Coated	100 yards	SP255 ● ●		SP257 ● ●			
	300 yards						SP2770 ● ●
Non-Absorbable Polyester, White Braided, Coated	100 yards			SP267 ● ●			

## LOOK™ Polyviolene™ Spools, Uncoated / Non-Sterile

Material	Length	Size						
		5-0	4-0	3-0	2-0	0	1	2
Non-Absorbable Polyester, Green Braided, Uncoated	100 yards	SP201 ● ●		SP203 ● ●	SP204 ● ●		SP206 ● ●	SP207 ● ●
Non-Absorbable Polyester, White Braided, Uncoated	100 yards	SP211 ● ●	SP210 ● ●					

Product numbers are color coded respectively as LOOK™ (purple), comparable Ethicon™ (black), and comparable Syneture™ (dark green).

● indicates that other manufacturers do not have a comparable product code or needle.

LOOK™ suture products are packaged in boxes of 12. LOOK™ spools and bulk suture are packaged as noted.

# Product Index

Code	Size	Material	Page Number
0F762	3-0	Silk, Black Braided Strands	27
0F766	3-0	Silk, Black Braided Strands	27
0F767	2-0	Silk, Black Braided Strands	27
0F770	2-0	Silk, Black Braided Strands	27
0F772	0	Silk, Black Braided Strands	27
0F782	3-0	Silk, Black Braided Strands	27
0F927	2-0	Nylon, Black Monofilament Strands	28
0F928	3-0	Nylon, Black Monofilament Strands	28
1011B	6-0	Polypropylene, Blue Monofilament	25
1013B	5-0	Polypropylene, Blue Monofilament	25
1018B	5-0	Polypropylene, Blue Monofilament	25
1020B	4-0	Polypropylene, Blue Monofilament	24
1022B	4-0	Polypropylene, Blue Monofilament	25
1031B	4-0	Polypropylene, Blue Monofilament	25
1058B	5-0	Polypropylene, Blue Monofilament	24
1075B	6-0	Polypropylene, Blue Monofilament	24
1077B	6-0	Polypropylene, Blue Monofilament	24
1078B	5-0	Polypropylene, Blue Monofilament	24
1079B	6-0	Polypropylene, Blue Monofilament	24
1082B	6-0	Polypropylene, Blue Monofilament	24
1085B	5-0	Polypropylene, Blue Monofilament	24
1091B	4-0	Polypropylene, Blue Monofilament	24
1095B	4-0	Polypropylene, Blue Monofilament	25
1176B	5-0	Polyviolene™ / Polyester, Green Braided Uncoated	26
1180B	4-0	Polyviolene™ / Polyester, Green Braided Uncoated	26
1205B	4-0	Polyviolene™ / Polyester, White Braided Uncoated	26
1241B	4-0	Plain Gut	9
1242B	6-0	Chromic Gut	11
1248B	5-0	Chromic Gut	11
1253B	4-0	Chromic Gut	11
1260B	6-0	Silk, Black Braided	18
1265B	5-0	Silk, Black Braided	18
1270B	4-0	Silk, Black Braided	18
1273B	6-0	Nylon, Black Monofilament	21
1275B	5-0	Nylon, Clear Monofilament	21
1276B	7-0	Nylon, Black Monofilament	21
1279B	6-0	Nylon, Black Monofilament	21
1280B	6-0	Nylon, Black Monofilament	21
1281B	6-0	Nylon, Black Monofilament	21
1282B	6-0	Nylon, Black Monofilament	21

# Product Index

Code	Size	Material	Page Number
1283B	5-0	Nylon, Black Monofilament	21
1284B	5-0	Nylon, Black Monofilament	21
1285B	5-0	Nylon, Black Monofilament	21
1288B	5-0	Nylon, Black Monofilament	21
1290B	4-0	Nylon, Black Monofilament	21
1291B	4-0	Nylon, Black Monofilament	21
1294B	3-0	Nylon, Black Monofilament	21
1415B	6-0	Plain Gut	9
1454B	4-0	Nylon, Black Monofilament	21
1456B	6-0	Nylon, Black Monofilament	21
310B	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
310BC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
311B	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
311BC	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
316B	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
316BC	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
317B	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
317BC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
332B	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
333B	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
339B	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
340B	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
341B	1	PolySyn™ / Polyglycolic Acid, Violet Braided	17
352B	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
353B	1	PolySyn™ / Polyglycolic Acid, Violet Braided	17
385B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
386B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
421B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
422B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
422BC	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
423B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
441B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
442B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
443B	2-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
480B	1	PolySyn™ / Polyglycolic Acid, Violet Braided	15
492B	6-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
493B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
494B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
495B	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
496B	4-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14

# Product Index

Code	Size	Material	Page Number
497B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
500B	5-0	Plain Gut	9
521B	5-0	Plain Gut	9
525B	3-0	Chromic Gut	11
531B	4-0	Plain Gut	9
535B	5-0	Plain Gut	9
536B	4-0	Plain Gut	9
537B	4-0	Plain Gut	10
543B	5-0	Chromic Gut	11
544B	4-0	Chromic Gut	12
547B	4-0	Chromic Gut	12
550B	5-0	Plain Gut	9
551B	4-0	Plain Gut	9
553B	3-0	Plain Gut	10
554B	3-0	Plain Gut	9
555B	2-0	Plain Gut	10
557B	5-0	Chromic Gut	11
558B	4-0	Chromic Gut	11
559B	4-0	Chromic Gut	11
560B	3-0	Chromic Gut	11
561B	2-0	Chromic Gut	11
562B	0	Chromic Gut	12
563B	5-0	Chromic Gut	11
564B	4-0	Chromic Gut	11
589B	4-0	Chromic Gut	12
590B	5-0	Chromic Gut	12
591B	4-0	Plain Gut	9
592B	3-0	Plain Gut	9
593B	3-0	Plain Gut	10
594B	3-0	Plain Gut	10
595B	4-0	Chromic Gut	11
596B	5-0	Plain Gut	9
598B	5-0	Chromic Gut	12
599B	4-0	Plain Gut	10
683B	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	14
750B	6-0	Silk, Black Braided	19
752B	6-0	Silk, Black Braided	18
754B	5-0	Silk, Black Braided	19
765B	3-0	Silk, Black Braided	20
768B	3-0	Silk, Black Braided	19

# Product Index

Code	Size	Material	Page Number
770B	2-0	Silk, Black Braided	20
772B	0	Silk, Black Braided	20
774B	5-0	Silk, Black Braided	18
775B	3-0	Silk, Black Braided	19
776B	4-0	Silk, Black Braided	19
777B	4-0	Silk, Black Braided	20
778B	3-0	Silk, Black Braided	20
779B	5-0	Silk, Black Braided	20
780B	4-0	Silk, Black Braided	20
781B	4-0	Silk, Black Braided	19
782B	3-0	Silk, Black Braided	20
783B	3-0	Silk, Black Braided	20
784B	3-0	Silk, Black Braided	19
785B	3-0	Silk, Black Braided	20
786B	3-0	Silk, Black Braided	19
787B	4-0	Silk, Black Braided	18
792B	4-0	Silk, Black Braided	18
795B	3-0	Silk, Black Braided	18
797B	5-0	Polyviolene™ / Polyester, Green Braided Uncoated	26
911B	6-0	Nylon, Black Monofilament	22
913B	5-0	Nylon, Black Monofilament	22
915B	6-0	Nylon, Black Monofilament	22
916B	6-0	Nylon, Black Monofilament	22
917B	5-0	Nylon, Black Monofilament	22
918B	5-0	Nylon, Black Monofilament	22
919B	5-0	Nylon, Black Monofilament	22
920B	4-0	Nylon, Black Monofilament	22
921B	4-0	Nylon, Black Monofilament	22
922B	4-0	Nylon, Black Monofilament	22
923B	3-0	Nylon, Black Monofilament	23
925B	3-0	Nylon, Black Monofilament	22
927B	2-0	Nylon, Black Monofilament	22
928B	3-0	Nylon, Black Monofilament	22
931B	4-0	Nylon, Black Monofilament	22
933B	2-0	Nylon, Black Monofilament	23
944B	4-0	Nylon, Black Monofilament	22
945B	5-0	Nylon, Black Monofilament	22
946B	5-0	Nylon, Clear Monofilament	22
947B	4-0	Nylon, Black Monofilament	22
950B	3-0	Nylon, Clear Monofilament	22

# Product Index

Code	Size	Material	Page Number
FA-893NC	5-0	Polyviolene™ / Polyester, White Braided Uncoated	26
GC-304NC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	16
GC-322BC	3-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-323NC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-354NC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-358NC	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-359NC	1	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-433NC	5-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	16
GC-451NC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	15
GC-452NC	3-0	PolySyn™ / Polyglycolic Acid, Undyed Braided	15
GC-518NC	0	PolySyn™ / Polyglycolic Acid, Violet Braided	17
GC-532NC	2-0	PolySyn™ / Polyglycolic Acid, Violet Braided	15
GC-936NC	0	PolySyn™ / Polyglycolic Acid, Violet Braided	16
GC-981NC	6-0	PolySyn™ / Polyglycolic Acid, Violet Braided	16
MC-397NC	4-0	PolySyn™ / Polyglycolic Acid, Violet Braided	14
MC-423NC	3-0	PolySyn FA™ / Polyglycolic Acid, Undyed Braided	13
MC-442NC	3-0	PolySyn FA™ / Polyglycolic Acid, Undyed Braided	13
MC-917NC	2-0	PolySyn FA™ / Polyglycolic Acid, Undyed Braided	13
SP102	6-0	Silk, Black Braided Spool, 25 yards	29
SP103	5-0	Silk, Black Braided Spool, 25 yards	29
SP103	4-0	Silk, Black Braided Spool, 25 yards	29
SP105	3-0	Silk, Black Braided Spool, 25 yards	29
SP106	2-0	Silk, Black Braided Spool, 25 yards	29
SP114	6-0	Silk, Black Braided Spool, 100 yards	29
SP115	5-0	Silk, Black Braided Spool, 100 yards	29
SP116	4-0	Silk, Black Braided Spool, 100 yards	29
SP117	3-0	Silk, Black Braided Spool, 100 yards	29
SP118	2-0	Silk, Black Braided Spool, 100 yards	29
SP119	0	Silk, Black Braided Spool, 100 yards	29
SP120	1	Silk, Black Braided Spool, 100 yards	29
SP121	2	Silk, Black Braided Spool, 50 yards	29
SP1230	3-0	Silk, Black Braided Spool, 2000 yards	29
SP135	5-0	Polypropylene, Blue Monofilament, 100 yards	29
SP136	4-0	Polypropylene, Blue Monofilament, 100 yards	29
SP137	3-0	Polypropylene, Blue Monofilament, 100 yards	29
SP138	2-0	Polypropylene, Blue Monofilament, 100 yards	29
SP175	5-0	Nylon, Black Braided Spool, 100 yards	28
SP176	4-0	Nylon, Black Braided Spool, 100 yards	28
SP185	5-0	Nylon, Black Monofilament Spool, 100 yards	28
SP186	4-0	Nylon, Black Monofilament Spool, 100 yards	28

# Product Index

Code	Size	Material	Page Number
SP187	3-0	Nylon, Black Monofilament Spool, 100 yards	28
SP188	2-0	Nylon, Black Monofilament Spool, 100 yards	28
SP189	0	Nylon, Black Monofilament Spool, 100 yards	28
SP190	1	Nylon, Black Monofilament Spool, 100 yards	28
SP201	5-0	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP203	3-0	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP204	2-0	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP206	1	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP207	2	Polyviolene™ / Polyester, Green Braided, Uncoated, Spool, 100 yards	30
SP210	4-0	Polyviolene™ / Polyester, White Braided, Uncoated, Spool, 100 yards	30
SP211	5-0	Polyviolene™ / Polyester, White Braided, Uncoated, Spool, 100 yards	30
SP252	4-0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 250 yards	30
SP2520	4-0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 3281 yards	30
SP252R	4-0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 250 yards	30
SP2555	0	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 100 yards	30
SP257	2	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 100 yards	30
SP267	2	Polyviolene™ / Polyester, White Braided, Coated, Spool, 100 yards	30
SP2770	5	Polyviolene™ / Polyester, Green Braided, Coated, Spool, 300 yards	30

Visit [SurgicalSpecialties.com](http://SurgicalSpecialties.com) And Click "Find A Product" For SKUs And Cross-Reference Information.

---

Surgical Specialties Corporation, 50 Braintree Hill Office Park, Suite 101, Braintree, MA 02184, U.S.A.  
Phone: 610.404.3520 • 877.991.1110 (USA and Canada only) • Fax: 610.404.4010  
All trademarks are property of their respective owners. ©2015 Surgical Specialties Corporation.  
All Rights Reserved. SSC-027 R0 5/15

**SURGICAL  
SPECIALTIES**  
CORPORATION™