

Choose Excellence,  
Expect Value

# SharpPoint<sup>™</sup> **PLUS**

Surgical Sutures

# Suture Catalog



# Choose Excellence, Expect Value

Surgical Specialties Corporation is the world's largest independent suture manufacturer. While headquartered in Westwood, Massachusetts, SSC is a global company with development and manufacturing facilities in the United States, Germany, Mexico, China, and the United Kingdom.

We design and manufacture high performance wound closure products including Quill®, Sharpoint™, and LOOK™ brand sutures.

We are pleased to offer the new Sharpoint™ PLUS suture portfolio featuring German engineered needles.

Our philosophy is that performance doesn't have to be expensive.



# Table of Contents

Suture Material Selection .....	3-4
Advanced Needle Technology .....	5
Needle Anatomy and Terminology .....	6
Needle Points and Cutting Edges .....	7
Needle Charts.....	8
<b>Absorbable Sutures</b>	
Express Gut.....	9
Plain Gut.....	9-10
Chromic Gut .....	11-12
PolySyn FA™ (polyglycolic acid) .....	13
Monoderm™ (glycolide and ε-caprolactone) .....	14-15
PolySyn™ (polyglycolic acid).....	16-19
PDO™ (polydioxanone) .....	20-21
<b>Non-Absorbable Sutures</b>	
Silk.....	22-23
Nylon™ .....	24-25
Polypropylene.....	26-29
Polyviolen™ XL (polyester)-Coated.....	30-31
<b>Cross Reference Guide</b>	
Ethicon Cross Reference.....	32-41
Covidien Cross Reference.....	42-46



# Guide to Suture Materials

## Absorbable

Material	Material Description	BSR Profile	Mass Absorption
<b>EXPRESS GUT™</b>  Comparable to: Fast Absorbing Gut Suture	Twisted <ul style="list-style-type: none"> <li>Natural</li> <li>Packaged in Wetting Solution</li> </ul> Common uses: Skin Closure	In vivo strength retention: 40% at 5 days 0% at 10 days  In-Vivo Study performed by NAMSA Data on file at Surgical Specialties Corporation	Essentially complete by 30 days  In-Vivo Study performed by NAMSA Data on file at Surgical Specialties Corporation
<b>GUT</b>  Comparable to: Surgical Gut Suture - Plain	Twisted <ul style="list-style-type: none"> <li>Plain &amp; Chromic</li> <li>Packaged in Wetting Solution</li> </ul> Common uses: General Closure, OB/Gyn, Urology, Bowel Anastomosis	In vivo strength retention: Varies	Varies depending upon type of suture (plain vs. chromic) and site of placement
<b>CHROMIC GUT</b>  Comparable to: Surgical Gut Suture - Chromic	Twisted <ul style="list-style-type: none"> <li>Plain &amp; Chromic</li> <li>Packaged in Wetting Solution</li> </ul> Common uses: General Closure, OB/Gyn, Urology, Bowel Anastomosis	In vivo strength retention: Varies	Varies depending upon type of suture (plain vs. chromic) and site of placement
<b>POLYSYN FA™</b> (polyglycolic acid)  · VICRYL RAPIDE™ * Velosorb™ Fast	Braided <ul style="list-style-type: none"> <li>Undyed (clear)</li> <li>Coated</li> </ul> Common uses: General Closure, OB/Gyn, Urology	In vivo strength retention: 48% at 7 days 0% at 14 days	Essentially complete between 49 and 63 days
<b>MONODERM™</b> (polyglycolic acid and e-caprolactone)  · MONOCRYL™ * Caprosyn™/Biosyn™	Monofilament <ul style="list-style-type: none"> <li>Undyed (clear)</li> <li>Dyed (violet)</li> </ul> Common uses: Skin Closure, OB/Gyn, Urology, Gastrointestinal	In vivo strength retention: 42% to 76% at 7 days 36% to 52% at 14 days	Essentially complete by 90 days
<b>POLYSYN™</b> (polyglycolic acid)  · VICRYL™ * Polysorb™/Dexon™	Braided <ul style="list-style-type: none"> <li>Undyed (clear)</li> <li>Dyed (violet)</li> </ul> Monofilament (9-0 & 10-0 only) <ul style="list-style-type: none"> <li>Undyed (clear)</li> <li>Dyed (violet)</li> <li>Coated &amp; Uncoated</li> </ul> Common uses: Skin Closure, OB/Gyn, Urology, Gastrointestinal	In vivo strength retention: 50% at 14 days 20% at 21 days	Essentially complete between 50 and 90 days
<b>PDO™</b> (polydioxanone)  · PDS II™ * Maxon™	Monofilament <ul style="list-style-type: none"> <li>Dyed (violet)</li> </ul> Common uses: Plastic Surgery, Orthopaedic, OB/Gyn, Gastrointestinal	In vivo strength retention: 80% to 90% at 14 days 60% to 82% at 28 days 47% to 79% at 42 days	Essentially complete between 180 and 220 days

· Comparable Ethicon™ Material \* Comparable Covidien™ Material


## Non-Absorbable

Material	Material Description	Qualities
<b>SILK</b>  · PERMA-HAND™ Silk * SofSilk™	<b>Braided</b> <ul style="list-style-type: none"> <li>Undyed (white)</li> <li>Dyed (black)</li> <li>Coated</li> </ul> <b>Common uses:</b> General & Skin Closure, Gastrointestinal, Cardiovascular Surgery, Plastic Surgery, Ophthalmology	<ul style="list-style-type: none"> <li>Excellent handling and tying characteristics</li> </ul>
<b>NYLON™</b>  · ETHILON™ * Monosof™/ Dermalon™	<b>Monofilament</b> <ul style="list-style-type: none"> <li>Undyed (clear)</li> <li>Dyed (blue and black)</li> <li>Packaged dry or in wetting solution</li> </ul> <b>Common uses:</b> General & Skin Closure, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery, Gastrointestinal	<ul style="list-style-type: none"> <li>Passes easily through tissue</li> </ul>
<b>POLYPROPYLENE</b>  · PROLENE™ * Surgipro™/Surgipro™ II	<b>Monofilament</b> <ul style="list-style-type: none"> <li>Dyed (blue)</li> </ul> <b>Common uses:</b> General & Skin Closure, Neurosurgery, Plastic Surgery, Ophthalmology, Microsurgery, Cardiovascular Surgery	<ul style="list-style-type: none"> <li>Inert</li> <li>Easy to handle</li> <li>Passes easily through tissue</li> </ul>
<b>POLYVIOLENE™ (UNCOATED)</b> <small>(poly (ethylene-terephthate))</small>  · MERSILENE™ * Surgidac™	<b>Braided</b> <ul style="list-style-type: none"> <li>Undyed (white)</li> <li>Dyed (green)</li> <li>Uncoated</li> </ul> <b>Common uses:</b> General & Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures	<ul style="list-style-type: none"> <li>Smooth</li> <li>Easy to handle</li> <li>Good knot security</li> </ul>
<b>POLYVIOLENE™ XL (COATED)</b>  · ETHIBOND™ * Tricon™	<b>Braided</b> <ul style="list-style-type: none"> <li>Undyed (white)</li> <li>Dyed (green)</li> <li>Coated</li> </ul> <b>Common uses:</b> General & Skin Closure, Plastic Surgery, Ophthalmology, Cardiovascular Surgery, Neurological Procedures	<ul style="list-style-type: none"> <li>Smooth</li> <li>Easy to handle</li> <li>Good knot security</li> </ul>

· Comparable Ethicon™ Material \* Comparable Covidien™ Material


# Advanced Needle Technology

Sharpoint Plus™ needle performance is an output of its design. Materials, geometry, and coatings work together to produce a high performance surgical instrument.

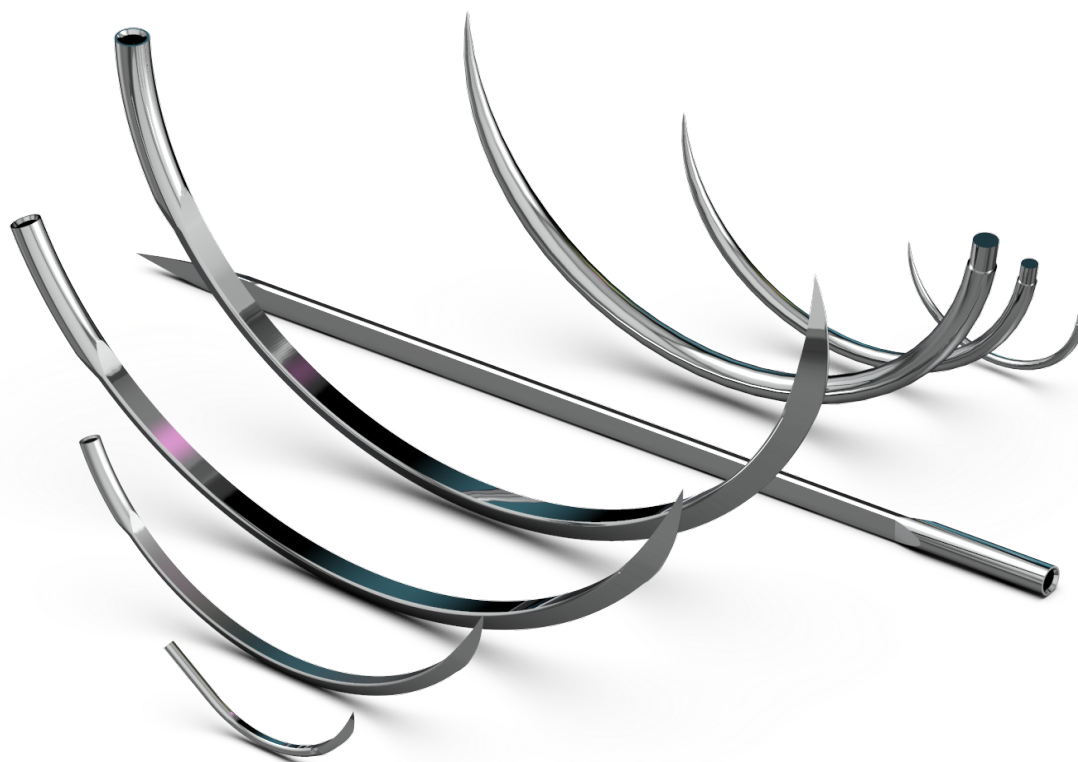


## Proprietary Stainless Steel Alloy


Sharpoint Plus™ needles are forged from an alloy of Iron, Nickel, and Chromium. Each element serves a role in delivering strength, ductility, and hardness.

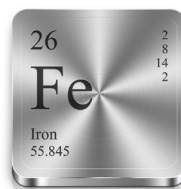

Strength

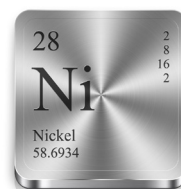

Ductility

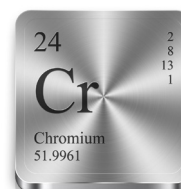

Hardness

# Needle Anatomy And Terminology

## Needle Terminology

**Length:** Distance measured along the circumference of the needle, from the point to the swage. It is depicted on the packaging.

**Shape:** Different surgeries or surgical techniques demand different needle shapes.

**Point:** The sharpened part of the needle that first penetrates the tissue. The point consists of the area from the sharpened tip of the needle to the widest cross section of the needle body. The point determines ease of penetration and the initial size and shape of the hole.


**Swage:** Connection point of the suture and the needle. Micro teeth grip the suture material and prevent it from slipping out.

**Wire Diameter:** Gauge or thickness of the wire originally used to form the needle body.


**Body:** The part between the point and the swage of the needle.

**Included Angle (Curvature):** The fraction of the circle included in the curve of the needle. This is depicted on the packaging.


## Needle Anatomy

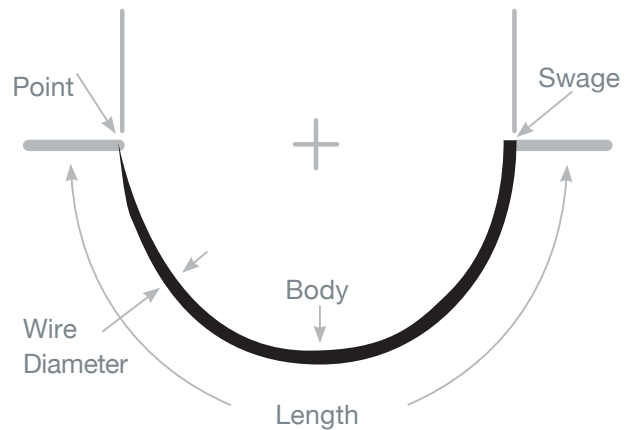


## Needle Shapes and General Uses

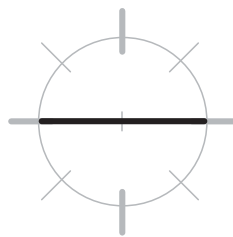


Straight (aka: Keith): Skin Closure

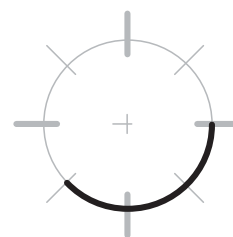

3/8 Circle: Skin Closure, Ophthalmic

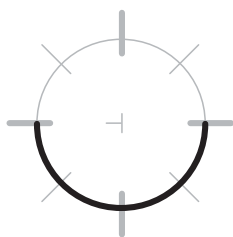



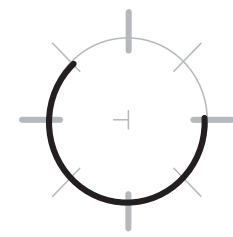


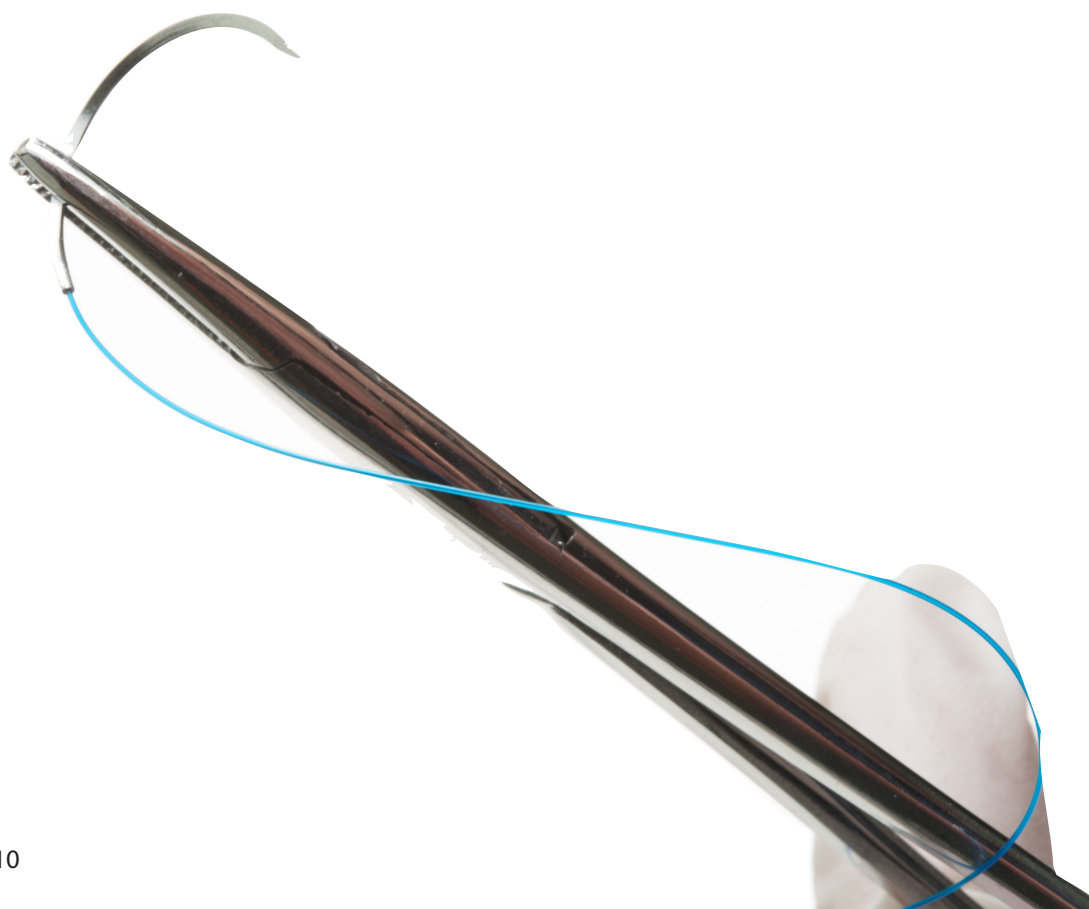
1/2 Circle: Skin Closure, GI, Cardiovascular, Ophthalmic



5/8 Circle: Robotics and Laparoscopic




# Needle Points And Cutting Edges

Needle Type	Point Geometry	Description	Commonly Used For
Precision Conventional Cutting		Triangular point with cutting edge on the inner curvature and precision-honed for added sharpness	General skin closure, subcutaneous tissue, sometimes for ophthalmic surgery, plastic or reconstructive surgery
Reverse Cutting		Triangular point with cutting edge on the outer curvature	Skin closure, ligamentous or fibrous tissues
Precision Reverse Cutting		Triangular point with cutting edge on the outer curvature and precision-honed for added sharpness	Primarily used in plastic surgery for delicate procedures and where a good cosmetic result is desired
Taper Point		Round shaft, straight or curved, taper point with no cutting edge	Soft tissue closure such as fascia, vascular, gastrointestinal, and most soft tissue below the skin
Taper Cutting		Four (4) cutting edges with a taper shaft	For penetration through tough tissue



























## Precision Conventional Cutting Needles

Surgical Specialties	Ethicon™	Covidien™	
3/8 Circle			
PC-1	PC-1	PC-13	
PC-3	PC-3	PC-11	
PC-5	PC-5	PC-12	



## Precision Reverse Cutting Needles

Surgical Specialties	Ethicon™	Covidien™	
3/8 Circle			
P-1	P-1	P-10	
G-1	G-1	G-1	
P-3	P-3	P-13	
G-3	G-3	G-3	
PS-3	PS-3	P-11	
PS-2	PS-2	P-12	
PS-1	PS-1	P-14	
PS	PS	P-15	
1/2 Circle			
PS-4	PS-4	P-24	








## Reverse Cutting Needles

Surgical Specialties	Ethicon™	Covidien™	
3/8 Circle			
C-2	C-2	C-2	
FS-3	FS-3	FS-3	
FS-2	FS-2	C-13	
FS-1	FS-1	C-14	
FS	FS	C-15	
FSL	FSL	C-16	
FSLX	FSLX	C-17	
			
1/2 Circle			
OS-4	OS-4	HOS-10	
CP-2	CP-2	GS-10	
OS-6	OS-6	HOS-11	
CP-1	CP-1	GS-11	
			





















## Straight Cutting Needles

Surgical Specialties	Ethicon™	Covidien™	
SC-1	SC-1		
KS	KS	SC-2	

## Taper Cutting Needles


Surgical Specialties	Ethicon™	Covidien™	
3/8 Circle			
CC	CC	KV-1	
CC-1	CC-1		
V-4	V-4	KV-4	
1/2 Circle			
V-5	V-5	KV-5	
V-7	V-7	KV-7	
V-34	V-34	KV-34	
			

## Taper Point Needles



Surgical Specialties	Ethicon™	Covidien™	
3/8 Circle			
BV-1	BV-1	CV-1	
BV	BV	CV	
C-1	C-1	CV-11	
BB	BB	CV-15	
1/2 Circle			
RB-2	RB-2	CV-22	
TF	TF	CVF-21	
RB-1	RB-1	CV-23	
SH-1	SH-1	CV-25	
SH	SH	V-20	
CT-2	CT-2	GS-22	
MO-6	MO-6	HGS-22	
MH	MH	V-26	
CT-1	CT-1	GS-21	
CT	CT	GS-24X	
CTX	CTX	GS-25	
			
			
5/8 Circle			
UR-6	UR-6	GU-46	
UR-5	UR-5	GU-45	
			

# Express Gut™

## Precision Conventional Cutting


				Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0	
<b>PC-1</b> 13mm 3/8 Circle 	Express Gut™	18"/45cm	B1916N	B1915N					

## Precision Reverse Cutting






				Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0	
<b>P-1</b> 11mm 3/8 Circle 	Express Gut™	18"/45cm	B1917N						
<b>P-3</b> 13mm 3/8 Circle 	Express Gut™	18"/45cm	B1918N	B1919N					

# Plain Gut

## Precision Conventional Cutting



				Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0	
<b>PC-1</b> 13mm 3/8 Circle 	Plain Gut	18"/45cm	B2711N						




## Precision Reverse Cutting




				Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0	
<b>G-1</b> 11mm 3/8 Circle 	Plain Gut	18"/45cm	B774N						
<b>G-1</b> 11mm 3/8 Circle 	Plain Gut	18"/45cm	B770N						
<b>P-3</b> 13mm 3/8 Circle 	Plain Gut	18"/45cm		B686N	B1644N				
<b>G-3</b> 13mm 3/8 Circle 	Plain Gut	18"/45cm		B772N					
<b>PS-2</b> 19mm 3/8 Circle 	Plain Gut	18"/45cm			B1627N	B592N			

N=box of 12; D=box of 36

# Plain Gut



Reverse Cutting 			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>FS-2</b> 19mm 3/8 Circle 	Plain Gut	27"/70cm			B821N	B822N		

Straight Cutting 			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>SC-1</b> 13mm Straight 	Plain Gut	18"/45cm			B1824N			
<b>SC-1</b> 13mm Straight 	Plain Gut	18"/45cm			B1828N			

Taper Point 			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>CT-1</b> 36mm 1/2 Circle 	Plain Gut	27"/30cm				B842N	B843N	
<b>CT</b> 40mm 1/2 Circle 	Plain Gut	27"/30cm				B852N	B853N	

# Chromic Gut








Precision Reverse Cutting 			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>G-1</b> 11mm 3/8 Circle 	Chromic Gut	12"/30cm	C796N					
<b>G-1</b> 11mm 3/8 Circle 	Chromic Gut	18"/45cm	C790N					
<b>P-3</b> 13mm 3/8 Circle 	Chromic Gut	18"/45cm		C687N	C1654N			
<b>G-3</b> 13mm 3/8 Circle 	Chromic Gut	18"/45cm		C792N	C793N			
<b>PS-3</b> 16mm 3/8 Circle 	Chromic Gut	18"/45cm		C1636N				
<b>PS-4</b> 16mm 1/2 Circle 	Chromic Gut	18"/45cm			C1643N			
<b>PS-2</b> 19mm 3/8 Circle 	Chromic Gut	18"/45cm			C1637N			
		27"/70cm				C1638N		

Reverse Cutting 			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>FS-2</b> 19mm 3/8 Circle 	Chromic Gut	27"/70cm			C635N	C636N		

N=box of 12; D=box of 36

11 surgicalspecialties.com | 877.991.1110

# Chromic Gut





Taper Point 			Size					
	Material	Length	5-0	4-0	3-0	2-0	0	1
<b>RB-1</b>  17mm 1/2 Circle	Chromic Gut	27"/70cm	C202N	C203N	C204N			
<b>SH-1</b>  22mm 1/2 Circle	Chromic Gut	27"/70cm		C181N	C182N			
<b>SH</b>  26mm 1/2 Circle	Chromic Gut	27"/70cm		C121N	C122N	C123N	C124N	
<b>CT-1</b>  36mm 1/2 Circle	Chromic Gut	27"/70cm			C810N	C811N	C812N	C813N
		36"/90cm			C922N	C923N	C924N	C925N
<b>CT</b>  40mm 1/2 Circle	Chromic Gut	27"/70cm					C802N	
		36"/90cm				C913N	C914N	C915N
<b>CTX</b>  48mm 1/2 Circle	Chromic Gut	36"/90cm					C904N	C905N

# PolySyn FA™ / Polyglycolic Acid

## Precision Reverse Cutting




### Size

	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>P-1</b>  11mm 3/8 Circle	Undyed Braided	18"/45cm		M490N				
<b>P-3</b>  13mm 3/8 Circle	Undyed Braided	18"/45cm		M493N	M494N			
<b>PS-2</b>  19mm 3/8 Circle	Undyed Braided	18"/45cm			M496N	M497N		
		27"/70cm			M426N			
<b>PS-1</b>  24mm 3/8 Circle	Undyed Braided	27"/70cm				M442N		

## Taper Cutting




### Size

	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>V-34</b>  36mm 1/2 Circle	Undyed Braided	36"/90cm					M917N	




## Taper Point






### Size

	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>CT-1</b>  36mm 1/2 Circle	Undyed Braided	36"/90cm				M944N	M945N	


Precision Conventional Cutting 

			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>PC-1</b> 13mm 3/8 Circle 	Undyed Monofilament	18"/45cm		Y834N				
<b>PC-3</b> 16mm 3/8 Circle 	Undyed Monofilament	18"/45cm		Y844N	Y845N			
<b>PC-5</b> 19mm 3/8 Circle 	Undyed Monofilament	18"/45cm			Y823N			

Precision Reverse Cutting 





			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>P-3</b> 13mm 3/8 Circle 	Undyed Monofilament	18"/45cm	Y492N	Y493N	Y494N			
	Violet Monofilament	18"/45cm		Y463N	Y464N			
<b>PS-2</b> 19mm 3/8 Circle 	Undyed Monofilament	18"/45cm		Y495N	Y496N	Y497N		
		27"/70cm			Y426N	Y427N		
	Violet Monofilament	18"/45cm			Y513N			
<b>PS-1</b> 24mm 3/8 Circle 	Undyed Monofilament	18"/45cm			Y682N			
		27"/70cm			Y935N	Y936N		

Reverse Cutting 

			Size				
	Material	Length	5-0	4-0	3-0	2-0	0
<b>CP-1</b> 36mm 1/2 Circle 	Undyed Monofilament	27"/70cm				Y266N	

# Monoderm™

Taper Point 




			Size					
	Material	Length	5-0	4-0	3-0	2-0	0	1
<b>RB-1</b>   17mm 1/2 Circle	Undyed Monofilament	27"/70cm	Y213D	Y214N	Y215N			
	Violet Monofilament	27"/70cm	Y303D	Y304D	Y305D			
<b>SH</b>   26mm 1/2 Circle	Undyed Monofilament	27"/70cm		Y415N	Y416N	Y417N Y417D		
	Violet Monofilament	27"/70cm			Y316N			
<b>CT-1</b>   36mm 1/2 Circle	Undyed Monofilament	36"/90cm			Y944N	Y945N	Y946N	
<b>CTX</b>   48mm 1/2 Circle	Violet Monofilament	36"/90cm					Y398N	Y399N

N=box of 12; D=box of 36





# PolySyn™/Polyglycolic Acid

## Precision Conventional Cutting

			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>PC-1</b> 13mm 3/8 Circle 	Undyed Braided	18"/45cm		G834N				
<b>PC-3</b> 16mm 3/8 Circle 	Undyed Braided	18"/45cm			G845N <i>ALT G835N</i>			
<b>PC-5</b> 19mm 3/8 Circle 	Undyed Braided	18"/45cm				G824N		






## Precision Reverse Cutting

			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>P-1</b> 11mm 3/8 Circle 	Undyed Braided	18"/45cm	G489N					
<b>P-3</b> 13mm 3/8 Circle 	Undyed Braided	18"/45cm	G492N	G493N	G494N			
	Violet Braided	18"/45cm		G463N	G464N			
<b>PS-3</b> 16mm 3/8 Circle 	Undyed Braided	18"/45cm		G500N				
<b>PS-4</b> 16mm 1/2 Circle 	Undyed Braided	18"/45cm		G506N	G507N			
<b>PS-2</b> 19mm 3/8 Circle 	Undyed Braided	18"/45cm		G495N	G496N	G497N		
		27"/70cm			G426N	G427N	G428N	
	Violet Braided	27"/70cm			G0392N			
<b>PS-1</b> 24mm 3/8 Circle 	Undyed Braided	18"/45cm				G683N		
		27"/70cm			G935N	G936N		

N=box of 12; D=box of 36

# PolySyn™/Polyglycolic Acid

Reverse Cutting ▾

			Size					
	Material	Length	5-0	4-0	3-0	2-0	0	1
<b>FS-2</b>  19mm 3/8 Circle	Undyed Braided	27"/70cm	G421N	G422N	G423N			
	Violet Braided	18"/45cm		G392N	G393N			
<b>FS-1</b>  24mm 3/8 Circle	Undyed Braided	27"/70cm		G441N	G442N	G443N		
	Violet Braided	27"/70cm		G451N				
<b>CP-2</b>  26mm 1/2 Circle	Undyed Braided	27"/70cm				G869N		
<b>CP-1</b>  36mm 1/2 Circle	Undyed Braided	27"/70cm				G266N	G267N	G268N
<b>OS-6</b>  36mm 1/2 Circle	Undyed Braided	27"/70cm						G535N


N=box of 12; D=box of 36

17surgicalspecialties.com | 877.991.1110

# PolySyn™/Polyglycolic Acid







## Taper Cutting



	Material	Length	Size				
			5-0	4-0	3-0	2-0	0
<b>V-34</b>  36mm 1/2 Circle	Undyed Braided	36"/90cm				G917N	






## Taper Point



	Material	Length	Size					
			5-0	4-0	3-0	2-0	0	1
<b>TF</b>  13mm 1/2 Circle	Undyed Braided	27"/70cm	G433N					
<b>RB-1</b>  17mm 1/2 Circle	Undyed Braided	27"/70cm		G214N	G215N			
<b>SH</b>  26mm 1/2 Circle	Undyed Braided	27"/70cm		G415N	G416N	G417N		
	Violet Braided	27"/70cm			G316N	G317N		
<b>SH-1</b>  22mm 1/2 Circle	Violet Braided	27"/70cm		G310N	G311N			
<b>CT-2</b>  26mm 1/2 Circle	Undyed Braided	27"/70cm				G269N	G270N	
	Violet Braided	27"/70cm				G333N	G334N	
<b>UR-6</b>  26mm 5/8 Circle	Violet Braided	27"/70cm				G602N	G603N	



N=box of 12; D=box of 36

# PolySyn™/Polyglycolic Acid



Taper Point 			Size					
			5-0	4-0	3-0	2-0	0	1
	Material	Length						
<b>CT-1</b>  36mm 1/2 Circle	Undyed Braided	27"/70cm			G258N	G259N	G260N	G261N
		36"/90cm			G944N	G945N	G946N	G947N
	Violet Braided	27"/70cm			G338N	G339N	G340N	G341N
		36"/90cm			G344N	G345N <i>ALT G354N</i>	G346N	G347N
<b>UR-5</b>  36mm 5/8 Circle	Violet Braided	27"/70cm					G376N	
<b>CT</b>  40mm 1/2 Circle	Undyed Braided	36"/90cm				G957N	G958N	
	Violet Braided	36"/90cm					G358N	
<b>CTX</b>  48mm 1/2 Circle	Undyed Braided	36"/90cm					G978N	G977N
	Violet Braided	36"/90cm					G370N	G371N

N=box of 12; D=box of 36




Precision Reverse Cutting 

			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>P-3</b> 13mm 3/8 Circle 	Violet Monofilament	18"/45cm		R493N	R494N			
<b>PS-2</b> 19mm 3/8 Circle 	Violet Monofilament	18"/45cm			R513N <i>ALT</i> <i>R496N</i>	R497N		







Reverse Cutting 










			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>FS-2</b> 19mm 3/8 Circle 	Violet Monofilament	27"/70cm			R422N	R423N		
<b>FS-1</b> 24mm 3/8 Circle 	Violet Monofilament	27"/70cm					R443N	






Taper Point 

			Size				
	Material	Length	7-0	6-0	5-0	4-0	3-0
<b>C-1</b> 13mm 3/8 Circle 	Violet Monofilament	30"/75cm		R127N	R126N		
<b>RB-2</b> 13mm 1/2 Circle 	Violet Monofilament	30"/75cm			R148N		
<b>RB-1</b> 17mm 1/2 Circle 	Violet Monofilament	27"/70cm			R303N	R304N	R305N

Taper Point 

			Size					
	Material	Length	5-0	4-0	3-0	2-0	0	1
<b>RB-1</b> 17mm 1/2 Circle 	Violet Monofilament	30"/75cm	R320N					
<b>SH</b> 26mm 1/2 Circle 	Violet Monofilament	27"/70cm			R316N	R317N		
<b>CT-2</b> 26mm 1/2 Circle 	Violet Monofilament	27"/70cm				R333N	R334N	
<b>CT-1</b> 36mm 1/2 Circle 	Violet Monofilament	27"/70cm				R339N	R340N	R341N
		36"/90cm					R346N	R347N
<b>CT</b> 40mm 1/2 Circle 	Violet Monofilament	36"/90cm					R358N	R359N
<b>CTX</b> 48mm 1/2 Circle 	Violet Monofilament	36"/90cm					R370N	R371N
		60"/120cm					R990N	


Precision Reverse Cutting 			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>P-1</b> 11mm 3/8 Circle 	Black Braided	18"/45cm	D639N					
<b>G-1</b> 11mm 3/8 Circle 	Black Braided	18"/45cm	D780N					
<b>P-3</b> 13mm 3/8 Circle 	Black Braided	18"/45cm	D1639N	D640N	D641N			
<b>G-3</b> 13mm 3/8 Circle 	Black Braided	18"/45cm			D783N			
<b>PS-1</b> 24mm 3/8 Circle 	Black Braided	18"/45cm				D1684N		
<b>PS-2</b> 19mm 3/8 Circle 	Black Braided	18"/45cm			D683N			
<b>PS</b> 26mm 3/8 Circle 	Black Braided	18"/45cm					D1588N	
<b>PSL</b> 30mm 3/8 Circle 	Black Braided	18"/45cm						D580N
		30"/75cm						D590N

Reverse Cutting 			Size					
	Material	Length	5-0	4-0	3-0	2-0	0	1
<b>FS-2</b> 19mm 3/8 Circle 	Black Braided	18"/45cm	D682N		D1679N			
		30"/75cm			D0768N			
<b>FS-1</b> 24mm 3/8 Circle 	Black Braided	18"/45cm			D684N	D0238N		
<b>FS</b> 26mm 3/8 Circle 	Black Braided	18"/45cm				D685N		
		30"/75cm			D0212N			
<b>FSL</b> 30mm 3/8 Circle 	Black Braided	18"/45cm					D678N	
		30"/75cm			D675N	D679N	D680N	




N=box of 12; D=box of 36

# Silk

## Straight Cutting

			Size					
No Needle	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>KS</b>  60mm Straight	Black Braided	30"/75cm				D622N	D623N D623D	

## Taper Point

			Size					
No Needle	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>C-1</b>  13mm 3/8 Circle	Black Braided	18"/45cm	D707N		D7734N			
<b>SH</b>  26mm 1/2 Circle	Black Braided	30"/75cm				D832N	D833N	D834N
<b>CT-1</b>  36mm 1/2 Circle	Black Braided	30"/75cm						D424N




## Silk Ligatures Tangle Free Dispenser

			Size					
No Needle	Material	Length	6-0	5-0	4-0	3-0	2-0	0
	Black Braided	12x18"/45cm			D183D	D184D	D185D	
		6x30"/75cm						D306D
		12x30"/75cm			D303D	D304D	D305D	







N=box of 12; D=box of 36



## Precision Conventional Cutting

			Size				
Needle	Material	Length	6-0	5-0	4-0	3-0	2-0
<b>PC-1</b> 13mm 3/8 Circle 	Black Monofilament	18"/45cm	A1856N	A1855N			
	Black Monofilament Pliable	18"/45cm		A1955N			
<b>PC-3</b> 16mm 3/8 Circle 	Black Monofilament	18"/45cm		A1865N	A1864N		
	Black Monofilament Pliable	18"/45cm		A1965N			
<b>PC-5</b> 19mm 3/8 Circle 	Black Monofilament	18"/45cm			A1894N	A1893N	
	Black Monofilament Pliable	18"/45cm			A1994N		




## Precision Reverse Cutting

			Size				
Needle	Material	Length	6-0	5-0	4-0	3-0	2-0
<b>P-1</b> 11mm 3/8 Circle 	Black Monofilament	18"/45cm	A697N				
<b>P-3</b> 13mm 3/8 Circle 	Black Monofilament	18"/45cm	A1698N	A698N <i>ALT AD-698N</i>	A699N		
	Undyed Monofilament	18"/45cm		A690N	A691N		
<b>PS-4</b> 16mm 1/2 Circle 	Black Monofilament	18"/45cm			A1662N		
<b>PS-3</b> 16mm 3/8 Circle 	Black Monofilament	18"/45cm	A1665N	A1668N			
<b>PS-2</b> 19mm 3/8 Circle 	Black Monofilament	18"/45cm	A1281N	A1666N <i>ALT AD-1666</i>	A1667N A662D	A1669N	A593N
	Blue Monofilament	18"/45cm			AG667N		
<b>PS-1</b> 24mm 3/8 Circle 	Black Monofilament	18"/45cm			A1629D	A1663N	








N=box of 12; D=box of 36

# Nylon™

## Precision Reverse Cutting

			Size				
Needle	Material	Length	6-0	5-0	4-0	3-0	2-0
<b>PS</b>  26mm 3/8 Circle	Black Monofilament	18"/45cm					A585N
<b>PSL</b>  30mm 3/8 Circle	Black Monofilament	30"/75cm				A1691N	
<b>PSLX</b>  36mm 3/8 Circle	Black Monofilament	30"/75cm				A1683N	A1697N

## Reverse Cutting

			Size				
	Material	Length	6-0	5-0	4-0	3-0	2-0
<b>C-2</b>  13mm 3/8 Circle	Black Monofilament	18"/45cm	A667N				
<b>FS-3</b>  16mm 3/8 Circle	Black Monofilament	18"/45cm	A0181N	A0556N			
<b>FS-2</b>  19mm 3/8 Circle	Black Monofilament	18"/45cm		A661N AC-0568D	A662N		
<b>FS-1</b>  24mm 3/8 Circle	Black Monofilament	18"/45cm			A1629N	A663N	
		30"/75cm				A669N AC-0586D	
<b>FS</b>  26mm 3/8 Circle	Black Monofilament	18"/45cm					A664N
<b>FSL</b>  30mm 3/8 Circle	Clear Monofilament	30"/75cm				A1671N	
<b>FSLX</b>  36mm 3/8 Circle	Black Monofilament	30"/75cm				A1673N	A1674N

N=box of 12; D=box of 36

Precision Conventional Cutting 				Size					
		Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>PC-1</b> 13mm 3/8 Circle		Blue Monofilament	18"/45cm	J8617N	J8618N				
<b>PC-3</b> 16mm 3/8 Circle		Blue Monofilament	18"/45cm		J8635N	J8634N			
<b>PC-5</b> 19mm 3/8 Circle		Blue Monofilament	18"/45cm			J8631N	J8632N		


Precision Reverse Cutting 				Size				
		Material	Length	7-0	6-0	5-0	4-0	3-0
<b>P-1</b> 11mm 3/8 Circle		Blue Monofilament	18"/45cm		J8697N			
<b>P-3</b> 13mm 3/8 Circle		Blue Monofilament	18"/45cm		J8695N	J8698N	J8699N	
<b>PS-3</b> 16mm 3/8 Circle		Blue Monofilament	18"/45cm		J8680N	J8681N		
<b>PS-2</b> 19mm 3/8 Circle		Blue Monofilament	18"/45cm			J8686N	J8682N	J8687N
<b>PS-1</b> 24mm 3/8 Circle		Blue Monofilament	18"/45cm					J8663N

Reverse Cutting 				Size				
		Material	Length	6-0	5-0	4-0	3-0	2-0
<b>FS-2</b> 19mm 3/8 Circle		Blue Monofilament	18"/45cm		J8661N	J8683N	J8665N	
<b>FS-1</b> 24mm 3/8 Circle		Blue Monofilament	18"/45cm				J8684N	
<b>FS</b> 26mm 3/8 Circle		Blue Monofilament	18"/45cm					J8685N





N=box of 12; D=box of 36

# Polypropylene


## Straight Cutting

			Size					
	Material	Length	6-0	5-0	4-0	3-0	2-0	0
<b>KS</b>  60mm Straight	Blue Monofilament	30"/75cm					J8623N	

## Taper Cutting







			Size					
	Material	Length	7-0	6-0	5-0	4-0	3-0	2-0
<b>CC</b>  9.3mm 3/8 Circle	Blue Monofilament	24"/60cm	J8704D	J8807D				
<b>CC-1</b>  13mm 3/8 Circle	Blue Monofilament	30"/75cm		J8707D				
<b>V-5</b>  17mm 1/2 Circle	Blue Monofilament	36"/90cm				J8935D		
<b>V-7</b>  26mm 1/2 Circle	Blue Monofilament	36"/90cm				J8975D	J8976D	J8977D

## Taper Point






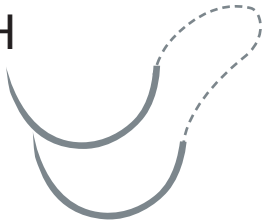



			Size					
	Material	Length	7-0	6-0	5-0	4-0	3-0	2-0
<b>BV-1</b>  9.3mm 3/8 Circle	Blue Monofilament	18"/45cm	J8701D	J8806D				
		24"/60cm	J8702D	J8805D				
			J8747N	J8807D				
			J8703D	J8709D				
		30"/75cm	J8755N					

N=box of 12; D=box of 36

## Taper Point

			Size					
	Material	Length	7-0	6-0	5-0	4-0	3-0	2-0
<b>C-1</b>  13mm 3/8 Circle	Blue Monofilament	18"/45cm		J8718D	J8717D			
		24"/60cm		J8726D	J8725D			
		30"/75cm		J8706D				
		36"/90cm			J8720D			
<b>BV</b>  11mm 3/8 Circle	Blue Monofilament	24"/60cm		J8610N				
		30"/75cm		J8776N				
<b>RB-2</b>  13mm 1/2 Circle	Blue Monofilament	18"/45cm		J8714D	J8713D			
		24"/60cm			J8554D			
		30"/75cm		J8711D	J8710D			
<b>BB</b>  17mm 3/8 Circle	Blue Monofilament	36"/90cm			J8580D	J8581D		
<b>RB-1</b>  17mm 1/2 Circle	Blue Monofilament	18"/45cm			J8756D			
		24"/60cm			J8555D			
		36"/90cm			J8556D	J8557D	J8558D	
<b>SH-1</b>  22mm 1/2 Circle	Blue Monofilament	30"/75cm				J8526D	J8762D	


# Polypropylene

Taper Point 		Size						
	Material	Length	5-0	4-0	3-0	2-0	0	1
<b>SH</b>  26mm 1/2 Circle	Blue Monofilament	30"/75cm			J8832N	J8833N	J8834N	
<b>SH</b>  26mm 1/2 Circle	Blue Monofilament	36"/90cm		J8521D	J8522D	J8523D		
		48"/122cm			J8534D	J8533D		
<b>MO-6</b>  26mm 1/2 Circle	Blue Monofilament	30"/75cm					J8418N	
<b>CT-2</b>  26mm 1/2 Circle	Blue Monofilament	30"/75cm				J8411N	J8412N	
<b>MH</b>  36mm 1/2 Circle	Blue Monofilament	36"/90cm			J8842D	J8843D		
		48"/122cm				J8853D		
<b>CT-1</b>  36mm 1/2 Circle	Blue Monofilament	30"/75cm				J8423N	J8424N	J8425N
<b>CT</b>  40mm 1/2 Circle	Blue Monofilament	30"/75cm					J8434N	J8435N
<b>CTX</b>  48mm 1/2 Circle	Blue Monofilament	30"/75cm						J8455N



N=box of 12; D=box of 36

# Polyviolene™ XL / Polyester - Coated





## Reverse Cutting ▾

		Size						
Material	Length	5-0	4-0	3-0	2-0	0	2	
<b>OS-4</b> 22mm 1/2 Circle 	Green Braided 30"/75cm						E519N	

## Taper Cutting ▾




Needle		Size						
Material	Length	6-0	5-0	4-0	3-0	2-0	0	
<b>V-4</b> 17mm 3/8 Circle 	Green Braided 30"/75cm						EX927D	
<b>V-7</b> 26mm 1/2 Circle 	Green Braided 30"/75cm						EX977D	

## Taper Point ⊙

		Size						
Material	Length	5-0	4-0	3-0	2-0	0	1	
<b>RB-1</b> 17mm 1/2 Circle 	Green Braided 30"/75cm		EX551D					
	36"/90cm	EX556D	EX557D	EX558D				
<b>SH-1</b> 22mm 1/2 Circle 	Green Braided 36"/90cm			EX762D	EX763D			
<b>SH</b> 26mm 1/2 Circle 	Green Braided 30"/75cm				EX563D			
	36"/90cm				EX523D	E524N		
	White Braided 36"/90cm				EX513D			
<b>SH</b> 26mm 1/2 Circle 	Green Braided 30"/75cm				E833N	E834N		

N=box of 12; D=box of 36

# Polyviolene™ XL / Polyester - Coated

Taper Point 		Material	Length	Size				
				5-0	4-0	3-0	2-0	0
<b>CT-2</b>  26mm 1/2 Circle	Green Braided	30"/75cm					E412N	
<b>CT-1</b>  36mm 1/2 Circle	Green Braided	30"/75cm					E424N	E425N

N=box of 12; D=box of 36



# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
424H	D424N	Silk 0 Black 30" CT-1 Taperpoint 1/2 Circle 36mm
580H	D580N	Silk 0 Black 18" PSL Reverse Cutting 3/8C 30mm
585H	A585N	Nylon™ 2-0 Black 18" PS Precision Reverse Cutting 3/8C 26mm
590H	D590N	Silk 0 Black 30" PSL Reverse Cutting 3/8C 30mm
593H	A593N	Nylon™ 2-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm
622H	D622N	Silk 3-0 Black 30" KS Reverse Cutting Straight 60mm
623H	D623D	Silk 2-0 Black 1x30" KS Straight Cutting 60mm
623H	D623N	Silk 2-0 Black 30" GS60 Straight Straight Cutting 60mm
635H	C635N	Chromic Gut 4-0 27" FS-2 Reverse Cutting 3/8C 19mm
636H	C636N	Chromic Gut 3-0 27" FS-2 Reverse Cutting 3/8C 19mm
639G	D639N	Silk 6-0 Black 18" P-1 Precision Reverse Cutting 3/8C 11mm
639G	D1639N	Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
640G	D640N	Silk 5-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm
641G	D641N	Silk 4-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm
661G	A661N	Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8C 19mm
661H	AC-0568D	Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8C 19mm
662G	A662N	Nylon™ 4-0 Black 18" FS-2 Reverse Cutting 3/8C 19mm
663G	A663N	Nylon™ 3-0 Black 18" FS-1 Reverse Cutting 3/8C 24mm
664G	A664N	Nylon™ 2-0 Black 18" FS Reverse Cutting 3/8C 26mm
667G	A667N	Nylon 6-0 Black 1x18" C-2 Reverse Cutting 3/8c 13mm
669H	A669N	Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8C 24mm
669H	AC-0586D	Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8C 24mm
675H	D675N	Silk 3-0 Black 1x30" FSL Reverse Cutting 3/8c 30mm
678G	D678N	Silk 0 Black 18" FSL Reverse Cutting 3/8C 30mm
679H	D679N	Silk 2-0 Black 30" FSL Reverse Cutting 3/8c 30mm
680H	D680N	Silk 0 Black 30" FSL Reverse Cutting 3/8C 30mm
682G	D682N	Silk, 5-0, Black, 18", FS-2, Reverse Cutting, 3/8 circle, 19mm
684G	D684N	Silk 3-0 Black 18" FS-1 Reverse Cutting 3/8C 24mm
685G	D685N	Silk 2-0 Black 18" FS Reverse Cutting 3/8C 26mm
686G	B686N	Plain Gut 5-0 18" P-3 Precision Reverse Cutting 3/8C 13mm
687G	C687N	Chromic Gut 5-0 18" P-3 Precision Reverse Cutting 3/8C 13mm
690G	A690N	Nylon™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
691G	A691N	Nylon™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
697H	A697N	Nylon™ 6-0 Black 18" P-1 Precision Reverse Cutting 3/8C 11mm
698G	AD-698N	Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm
698H	A698N	Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm
699G	A699N	Nylon™ 4-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm
707G	D707N	Silk 6-0 Black 18" C-1 Taper Point 13mm 3/8c
770G	B770N	Plain Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm Double Armed

# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
772G	B772N	Plain Gut 5-0 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed
774G	B774N	Plain Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm
780G	D780N	Silk 6-0 Black 18" G-1 Precision Reverse Cutting 3/8C 11mm Double Armed
783G	D783N	Silk 4-0 Black 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed
786G	D1639N	Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
790G	C790N	Chromic Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm Double Armed
792G	C792N	Chromic Gut 5-0 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed
793G	C793N	Chromic Gut 4-0 18" G-3 Precision Reverse Cutting 3/8C 13mm Double Armed
796G	C796N	Chromic Gut 6-0 18" G-1 Precision Reverse Cutting 3/8C 11mm
802H	C802N	Chromic Gut 0 27" CT Taperpoint 1/2 Circle 40mm
810H	C810N	Chromic Gut 3-0 27" CT-1 Taperpoint 1/2 Circle 36mm
811H	C811N	Chromic Gut 2-0 27" CT-1 Taperpoint 1/2 Circle 36mm
812H	C812N	Chromic Gut 0 27" CT-1 Taperpoint 1/2 Circle 36mm
813H	C813N	Chromic Gut 1 27" CT-1 Taperpoint 1/2 Circle 36mm
842H	B842N	Plain Gut 3-0 27" CT-1 Taperpoint 1/2 Circle 36mm
843H	B843N	Plain Gut 2-0 27" CT-1 Taperpoint 1/2 Circle 36mm
852H	B852N	Plain Gut 3-0 27" CT Taperpoint 1/2 Circle 40mm
853H	B853N	Plain Gut 2-0 27" CT Taperpoint 1/2 Circle 40mm
904H	C904N	Chromic Gut 0 36" CTX Taperpoint 1/2 Circle 48mm
905H	C905N	Chromic Gut 1 36" CTX Taperpoint 1/2 Circle 48mm
913H	C913N	Chromic Gut 2-0 36" CT Taperpoint 1/2 Circle 40mm
914H	C914N	Chromic Gut 0 36" CT Taperpoint 1/2 Circle 40mm
915H	C915N	Chromic Gut 1 36" CT Taperpoint 1/2 Circle 40mm
922H	C922N	Chromic Gut 3-0 36" CT-1 Taperpoint 1/2 Circle 36mm
923H	C923N	Chromic Gut 2-0 36" CT-1 Taperpoint 1/2 Circle 36mm
924H	C924N	Chromic Gut 0 36" CT-1 Taperpoint 1/2 Circle 36mm
925H	C925N	Chromic Gut 1 36" CT-1 Taperpoint 1/2 Circle 36mm
1588H	D1588N	Silk 2-0 Black 18" PS Precision Reverse Cutting 3/8C 26mm
1627H	B1627N	Plain Gut 4-0 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1629H	A1629N	Nylon™ 4-0 Black 18" FS-1 Reverse Cutting 3/8C 24mm
1636G	C1636N	Chromic Gut 5-0 18" PS-3 Precision Reverse Cutting 3/8C 16mm
1637G	C1637N	Chromic Gut 4-0 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1638H	C1638N	Chromic Gut 3-0 27" PS-2 Precision Reverse Cutting 3/8C 19mm
1639G	D1639N	Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
1643G	C1643N	Chromic Gut 4-0 18" PS-4 Precision Reverse Cutting 1/2C 16mm
1644G	B1644N	Plain Gut 4-0 18" P-3 Precision Reverse Cutting 3/8C 13mm
1654G	C1654N	Chromic Gut 4-0 18" P-3 Precision Reverse Cutting 3/8C 13mm
1662G	A1662N	Nylon™ 4-0 Black 18" PS-4 Precision Reverse Cutting 1/2C 16mm
1663G	A1663N	Nylon™ 3-0 Black 18" PS-1 Precision Reverse Cutting 3/8C 24mm
1665G	A1665N	Nylon™ 6-0 Black 18" PS-3 Precision Reverse Cutting 3/8C 16mm

# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
1666G	A1666N	Nylon™ 5-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1666G	AD-1666	Nylon™ 5-0 Black 1x18" PS-2 Precision Reverse Cutting 3/8c 19mm
1667G	A1667N	Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1667H	A662D	Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1668G	A1668N	Nylon™ 5-0 Black 18" PS-3 Precision Reverse Cutting 3/8C 16mm
1668G	A0556N	Nylon™ 5-0 Black 1x18" FS-3 Reverse Cutting 3/8c 16mm
1668H	A0556N	Nylon™ 5-0 Black 1x18" FS-3 Reverse Cutting 3/8c 16mm
1669H	A1669N	Nylon™ 3-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1671H	A1671N	Nylon™ 3-0 Black 30" FSL Reverse Cutting 3/8C 30mm
1673H	A1673N	Nylon™ 3-0 Black 30" FSLX Reverse Cutting 3/8C 36mm
1674H	A1674N	Nylon™ 2-0 Black 30" FSLX Reverse Cutting 3/8C 36mm
1677G	D683N	Silk 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8C 19mm
1679H	D1679N	Silk 3-0 Black 1x18" FS-2 Reverse Cutting 3/8c 19mm
1683H	A1683N	Nylon™ 3-0 Black 30" PSLX Reverse Cutting 3/8C 36mm
1690H	D675N	Silk 3-0 Black 30" FSL Precision Reverse Cutting 3/8C 30mm
1691H	A1691N	Nylon™ 3-0 Black 30" PSL Reverse Cutting 3/8C 30mm
1697H	A1697N	Nylon™ 2-0 Black 30" PSLX Reverse Cutting 3/8C 36mm
1698G	A1698N	Nylon™ 6-0 Black 18" P-3 Precision Reverse Cutting 3/8C 13mm
1824H	B1824N	Plain Gut 4-0 18" SC-1 Reverse Cutting Straight 13mm
1828H	B1828N	Plain Gut 4-0 18" SC-1 Reverse Cutting Straight 13mm Double Armed
1855G	A1855N	Nylon™ 5-0 Black 18" PC-1 Precision Conventional Cutting 3/8C 13mm
1856G	A1856N	Nylon™ 6-0 Black 18" PC-1 Precision Conventional Cutting 3/8C 13mm
1864G	A1864N	Nylon™ 4-0 Black 18" PC-3 Precision Conventional Cutting 3/8C 16mm
1865G	A1865N	Nylon™ 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8C 16mm
1893G	A1893N	Nylon™ 3-0 Black 18" PC-5 Precision Conventional Cutting 3/8C 19mm
1894G	A1894N	Nylon™ 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8C 19mm
1915G	B1915N	Express Gut™ 5-0 18" PC-1 Precision Conventional Cutting 3/8C 13mm
1916G	B1916N	Express Gut™ 6-0 18" PC-1 Precision Conventional Cutting 3/8C 13mm
1955G	A1955N	Nylon™ (pliable) 5-0 Black 18" PC-1 Precision Conventional Cutting 3/8C 13mm
1965G	A1965N	Nylon™ (pliable) 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8C 16mm
1994G	A1994N	Nylon™ (Pliable) 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8C 19mm
7734G	D7734N	Silk 4-0 Black 18" C-1 Taper Point 3/8c 13mm
8410H	J8832N	Polypropylene 3-0 Blue 30" SH Taperpoint 1/2 Circle 26mm
8411H	J8833N	Polypropylene 2-0 Blue 30" SH Taperpoint 1/2 Circle 26mm
8423H	J8423N	Polypropylene 2-0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm
8424H	J8424N	Polypropylene 0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm
8425H	J8425N	Polypropylene 1 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm
8434H	J8434N	Polypropylene 0 Blue 30" CT Taperpoint 1/2 Circle 40mm
8435H	J8435N	Polypropylene 1 Blue 30" CT Taperpoint 1/2 Circle 40mm
8455H	J8455N	Polypropylene 1 Blue 30" CTX Taperpoint 1/2 Circle 48mm

# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
8521H	J8521D	Polypropylene 4-0 Blue 36" SH Taperpoint 1/2C 26mm Double Armed
8522H	J8522D	Polypropylene 3-0 Blue 36" SH Taperpoint 1/2C 26mm Double Armed
8523H	J8523D	Polypropylene 2-0 Blue 36" SH Taperpoint 1/2C 26mm Double Armed
8526H	J8526D	Polypropylene 4-0 Blue 30" SH-1 Taperpoint 1/2C 22mm Double Armed
8534H	J8534D	Polypropylene 3-0 Blue 48" SH Taperpoint 1/2C 26mm Double Armed
8554H	J8554D	Polypropylene 5-0 Blue 24" RB-2 Taperpoint 1/2C 13mm Double Armed
8555H	J8555D	Polypropylene 5-0 Blue 24" RB-1 Taperpoint 1/2C 17mm Double Armed
8556H	J8556D	Polypropylene 5-0 Blue 36" RB-1 Taperpoint 1/2C 17mm Double Armed
8557H	J8557D	Polypropylene 4-0 Blue 36" RB-1 Taperpoint 1/2C 17mm Double Armed
8558H	J8558D	Polypropylene 3-0 Blue 36" RB-1 Taperpoint 1/2C 17mm Double Armed
8580H	J8580D	Polypropylene 5-0 Blue 36" BB Taperpoint 3/8C 17mm Double Armed
8581H	J8581D	Polypropylene 4-0 Blue 36" BB Taperpoint 3/8C 17mm Double Armed
8617G	J8617N	Polypropylene 6-0 Blue 18" PC-1 Precision Conventional Cutting 3/8C 13mm
8618G	J8618N	Polypropylene 5-0 Blue 18" PC-1 Precision Conventional Cutting 3/8C 13mm
8623H	J8623N	Polypropylene 2-0 Blue 1x30" KS Straight Cutting 60mm
8631G	J8631N	Polypropylene 4-0 Blue 18" PC-5 Precision Conventional Cutting 3/8C 19mm
8632G	J8632N	Polypropylene 3-0 Blue 18" PC-5 Precision Conventional Cutting 3/8C 19mm
8634G	J8634N	Polypropylene 4-0 Blue 18" PC-3 Precision Conventional Cutting 3/8C 16mm
8635G	J8635N	Polypropylene 5-0 Blue 18" PC-3 Precision Conventional Cutting 3/8C 16mm
8661G	J8661N	Polypropylene 5-0 Blue 18" FS-2 Reverse Cutting 3/8C 19mm
8663G	J8663N	Polypropylene 3-0 Blue 18" PS-1 Precision Reverse Cutting 3/8C 24mm
8665G	J8665N	Polypropylene 3-0 Blue 18" FS-2 Reverse Cutting 3/8C 19mm
8680G	J8680N	Polypropylene 6-0 Blue 18" PS-3 Precision Reverse Cutting 3/8C 16mm
8681G	J8681N	Polypropylene 5-0 Blue 18" PS-3 Precision Reverse Cutting 3/8C 16mm
8682G	J8682N	Polypropylene 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm
8683G	J8683N	Polypropylene 4-0 Blue 18" FS-2 Reverse Cutting 3/8C 19mm
8684G	J8684N	Polypropylene 3-0 Blue 18" FS-1 Reverse Cutting 3/8C 24mm
8685H	J8685N	Polypropylene 2-0 Blue 18" FS Reverse Cutting 3/8C 26mm
8686G	J8686N	Polypropylene 5-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm
8687H	J8687N	Polypropylene 3-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm
8695G	J8695N	Polypropylene 6-0 Blue 18" P-3 Precision Reverse Cutting 3/8C 13mm
8697G	J8697N	Polypropylene 6-0 Blue 18" P-1 Precision Reverse Cutting 3/8C 11mm
8698G	J8698N	Polypropylene 5-0 Blue 18" P-3 Precision Reverse Cutting 3/8C 13mm
8699G	J8699N	Polypropylene 4-0 Blue 18" P-3 Precision Reverse Cutting 3/8C 13mm
8701H	J8701D	Polypropylene 7-0 Blue 18" BV-1 Taperpoint 3/8C 9.3mm Double Armed
8702H	J8702D	Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8C 9.3mm Double Armed
8702H	J8747N	Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
8703H	J8703D	Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8C 9.3mm Double Armed
8704H	J8704D	Polypropylene 7-0 Blue 24" CC Tapercut 3/8C 9.3mm Double Armed
8706H	J8706D	Polypropylene 6-0 Blue 30" C-1 Taperpoint 3/8C 13mm Double Armed

# Ethicon Cross Reference Guide

Ethicon	Sharp point Plus™	Product Description
8707H	J8707D	Polypropylene 6-0 Blue 30" CC-1 Tapercut Cardio 3/8C 13mm Double Armed
8709H	J8709D	Polypropylene 6-0 Blue 30" BV-1 Taperpoint 3/8C 9.3mm Double Armed
8710H	J8710D	Polypropylene 5-0 Blue 30" RB-2 Taperpoint 1/2C 13mm Double Armed
8711H	J8711D	Polypropylene 6-0 Blue 30" RB-2 Taperpoint 1/2C 13mm Double Armed
8713H	J8713D	Polypropylene 5-0 Blue 18" RB-2 Taperpoint 1/2C 13mm Double Armed
8714H	J8714D	Polypropylene 6-0 Blue 18" RB-2 Taperpoint 1/2C 13mm Double Armed
8717H	J8717D	Polypropylene 5-0 Blue 18" C-1 Taperpoint 3/8C 13mm Double Armed
8718H	J8718D	Polypropylene 6-0 Blue 18" C-1 Taperpoint 3/8C 13mm Double Armed
8720H	J8720D	Polypropylene 5-0 Blue 36" C-1 Taperpoint 3/8C 13mm Double Armed
8725H	J8725D	Polypropylene 5-0 Blue 24" C-1 Taperpoint 3/8C 13mm Double Armed
8726H	J8726D	Polypropylene 6-0 Blue 24" C-1 Taperpoint 3/8C 13mm Double Armed
8755H	J8755N	Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
8756H	J8756D	Polypropylene 5-0 Blue 18" RB-1 Taperpoint 1/2C 17mm Double Armed
8762H	J8762D	Polypropylene 3-0 Blue 30" SH-1 Taperpoint 1/2C 22mm Double Armed
8805H	J8805D	Polypropylene 6-0 Blue 24" BV-1 Taperpoint 3/8C 9.3mm Double Armed
8806H	J8806D	Polypropylene 6-0 Blue 18" BV-1 Taperpoint 3/8C 9.3mm Double Armed
8807H	J8807D	Polypropylene 6-0 Blue 24" CC Tapercut 3/8C 9.3mm Double Armed
8833H	J8411N	Polypropylene 2-0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm
8834H	J8412N	Polypropylene 0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm
8834H	J8418N	Polypropylene 0 Blue 30" MO-6 Taperpoint 1/2 Circle 26mm
8834H	J8834N	Polypropylene 0 Blue 30" SH Taperpoint 1/2 Circle 26mm
8842H	J8842D	Polypropylene 3-0 Blue 36" MH Taperpoint 1/2C 36mm Double Armed
8843H	J8843D	Polypropylene 2-0 Blue 36" MH Taperpoint 1/2C 36mm Double Armed
8853H	J8853D	Polypropylene 2-0 Blue 48" MH Taperpoint 1/2C 36mm Double Armed
8935H	J8935D	Polypropylene 4-0 Blue 36" V-5 Tapercut Cardio 1/2C 17mm Double Armed
8975H	J8975D	Polypropylene 4-0 Blue 36" V-7 Tapercut Cardio 1/2C 26mm Double Armed
8976H	J8976D	Polypropylene 3-0 Blue 36" V-7 Tapercut Cardio 1/2C 26mm Double Armed
8977H	J8977D	Polypropylene 2-0 Blue 36" V-7 Tapercut Cardio 1/2C 26mm Double Armed
8886291451	EX513D	Polyviolene™ XL 2-0 36" SH Taper Point 1/2 Circle 26mm
8886295351	EX513D	Polyviolene™ XL 2-0 36" SH Taper Point 1/2 Circle 26mm
A183H	D183D	Silk 4-0 Black 18" Silk Ligatures No Needle
A184H	D184D	Silk 3-0 Black 18" Silk Ligatures No Needle
A185H	D185D	Silk 2-0 Black 18" Silk Ligatures No Needle
A303H	D303D	Silk 4-0 Black 30" Silk Ligatures No Needle
A304H	D304D	Silk 3-0 Black 30" Silk Ligatures No Needle
A305H	D305D	Silk 2-0 Black 30" Silk Ligatures No Needle
A306H	D306D	Silk 0 Black 30" Silk Ligatures No Needle
EP8610H	J8610N	Polypropylene 6-0 Blue 24" BV Taperpoint 3/8 Circle 11mm Double Armed
EP8776H	J8776N	Polypropylene 6-0 Blue 30" BV Taperpoint 3/8 Circle 11mm Double Armed
G121H	C121N	Chromic Gut 4-0 27" SH Taperpoint 1/2 Circle 26mm

# Ethicon Cross Reference Guide

Ethicon	Sharpoint Plus™	Product Description
G122H	C122N	Chromic Gut 3-0 27" SH Taperpoint 1/2 Circle 26mm
G123H	C123N	Chromic Gut 2-0 27" SH Taperpoint 1/2 Circle 26mm
G124H	C124N	Chromic Gut 0 27" SH Taperpoint 1/2 Circle 26mm
G181H	C181N	Chromic Gut 4-0 27" SH-1 Taperpoint 1/2 Circle 22mm
G182H	C182N	Chromic Gut 3-0 30" HR22 1/2 Circle Taper Point 22mm
G667G	AG667N	Nylon™ 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8C 19mm
H821H	B821N	Plain Gut 4-0 27" FS-2 Reverse Cutting 3/8C 19mm
H822H	B822N	Plain Gut 3-0 27" FS-2 Reverse Cutting 3/8C 19mm
J214H	G214N	PolySyn™ 4-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm
J215H	G215N	PolySyn™ 3-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm
J232H	G416N	PolySyn™ 3-0 Undyed 27" HR26 1/2 Circle Taper Point 26mm
J258H	G258N	PolySyn™ 3-0 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm
J259H	G259N	PolySyn™ 2-0 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm
J260H	G260N	PolySyn™ 0 Undyed 27" HR36 1/2 Circle Taper Point 36mm
J261H	G261N	PolySyn™ 1 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm
J266H	G266N	PolySyn™ 2-0 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm
J267H	G267N	PolySyn™ 0 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm
J268H	G268N	PolySyn™ 1 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm
J269H	G269N	PolySyn™ 2-0 Undyed 27" HR26 1/2 Circle Taper Point 26mm
J269H	G417N	PolySyn™ 2-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm
J270H	G270N	PolySyn™ 0 Undyed 27" CT-2 Taperpoint 1/2 Circle 26mm
J310H	G310N	Polysyn™ 4-0 Violet 27" SH-1 Taper Point 1/2c 22mm
J311H	G311N	Polysyn™ 3-0 Violet 1x27" SH-1 Taper Point 1/2c 22mm
J316H	G316N	PolySyn™ 3-0 Violet 27" HR26 1/2 Circle Taper Point 26mm
J317H	G317N	PolySyn™ 2-0 Violet 27" HR26 1/2 Circle Taper Point 26mm
J333H	G333N	PolySyn™ 2-0 Violet 27" HR26 1/2 Circle Taper Point 26mm
J334H	G334N	PolySyn™ 0 Violet 27" CT-2 Taperpoint 1/2 Circle 26mm
J338H	G338N	PolySyn™ 3-0 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm
J339H	G339N	PolySyn™ 2-0 Violet 27" HR36 1/2 Circle Taper Point 36mm
J340H	G340N	PolySyn™ 0 Violet 27" HRX36 1/2 Circle Taper Point 36.6mm
J341H	G341N	PolySyn™ 1 Violet 27" HRX36 1/2 Circle Taper Point 36.6mm
J344H	G344N	PolySyn™ 3-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
J345H	G345N	PolySyn™ 2-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
J345H	G354N	Polysyn™ 2-0 Violet 36" CT-1 Taper Point 1/2c 36mm
J346H	G346N	PolySyn™ 0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
J347H	G347N	PolySyn™ 1 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
J358H	G358N	PolySyn™ 0 Violet 36" HRX40 1/2 Circle Taper Point 40mm
J370H	G370N	PolySyn™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm
J371H	G371N	PolySyn™1 Violet 36" CTX Taperpoint 1/2 Circle 48mm
J376H	G376N	PolySyn™ 0 Violet 27" UR-5 Taperpoint 5/8C 36mm

# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
J392H	G392N	PolySyn™ 4-0 Violet 18" FS-2 Reverse Cutting 3/8C 19mm
J393H	G393N	PolySyn™ 3-0 Violet 18" FS-2 Reverse Cutting 3/8C 19mm
J415H	G415N	PolySyn™ 4-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm
J421H	G421N	Polysyn™ 5-0 Undyed 27" FS-2 Reverse Cutting 3/8c 19mm
J422H	G422N	PolySyn™ 4-0 Undyed 27" FS-2 Reverse Cutting 3/8C 19mm
J423H	G423N	PolySyn™ 3-0 Undyed 27" FS-2 Reverse Cutting 3/8C 19mm
J426H	G426N	PolySyn™ 4-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm
J427H	G427N	PolySyn™ 3-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm
J428H	G428N	PolySyn™ 2-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm
J433H	G433N	Polysyn™ 5-0 Undyed 27" TF Taper Point 1/2c 13mm
J441H	G441N	PolySyn™ 4-0 Undyed 27" FS-1 Reverse Cutting 3/8C 24mm
J442H	G442N	PolySyn™ 3-0 Undyed 27" FS-1 Reverse Cutting 3/8C 24mm
J443H	G443N	PolySyn™ 2-0 Undyed 27" FS-1 Reverse Cutting 3/8C 24mm
J451H	G451N	Polysyn™ 4-0 Violet 27" PS-1 Precision Reverse Cutting 3/8c 24mm
J463G	G463N	PolySyn™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm
J464G	G464N	PolySyn™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm
J489G	G489N	PolySyn™ 6-0 Undyed 18" P-1 Precision Reverse Cutting 3/8C 11mm
J492G	G492N	PolySyn™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
J493G	G493N	PolySyn™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
J494G	G494N	PolySyn™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
J495G	G495N	PolySyn™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
J496G	G496N	PolySyn™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
J497G	G497N	PolySyn™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
J500G	G500N	PolySyn™ 5-0 Undyed 18" PS-3 Precision Reverse Cutting 3/8C 16mm
J506G	G506N	Polysyn, 5-0, Undyed, 18", PS-4, Precision Reverse Cutting, 1/2 Circle, 16mm
J507G	G507N	PolySyn™ 4-0 Undyed 18" PS-4 Precision Reverse Cutting 1/2C 16mm
J535H	G535N	PolySyn™ 1 Undyed 27" OS-6 Reverse Cutting 1/2C 36mm
J602H	G602N	PolySyn™ 2-0 Violet 27" UR-6 Taperpoint 5/8C 26mm
J603H	G603N	PolySyn™ 0 Violet 27" UR-6 Taperpoint 5/8C 26mm
J683G	G683N	PolySyn™ 3-0 Undyed 18" PS-1 Precision Reverse Cutting 3/8C 24mm
J824G	G824N	PolySyn™ 3-0 Undyed 18" PC-5 Precision Conventional Cutting 3/8C 19mm
J834G	G834N	PolySyn™ 5-0 Undyed 18" PC-1 Precision Conventional Cutting 3/8C 13mm
J845G	G835N	Polysyn™ 4-0 Undyed 18" Precision Conventional Cutting 3/8 16mm
J845G	G845N	PolySyn™ 4-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8C 16mm
J869H	G869N	PolySyn™ 2-0 Undyed 27" CP-2 Reverse Cutting 1/2C 26mm
J917H	G917N	PolySyn™ 2-0 Undyed 36" V-34 Tapercut 1/2 Circle 36mm
J935H	G935N	PolySyn™ 4-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm
J936H	G936N	PolySyn™ 3-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm
J944H	G944N	PolySyn™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
J945H	G945N	PolySyn™ 2-0 Undyed 36" HR36 1/2 Circle Taper Point 36mm

# Ethicon Cross Reference Guide

Ethicon	Sharpoint Plus™	Product Description
J946H	G946N	PolySyn™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
J947H	G947N	PolySyn™ 1 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
J957H	G957N	PolySyn™ 2-0 Undyed 36" CT Taperpoint 1/2 Circle 40mm
J958H	G958N	PolySyn™ 0 Undyed 36" CT Taperpoint 1/2 Circle 40mm
J977H	G977N	PolySyn™ 1 Undyed 36" CTX Taperpoint 1/2 Circle 48mm
J978H	G978N	PolySyn™ 0 Undyed 36" CTX Taperpoint 1/2 Circle 48mm
K832H	D832N	Silk 3-0 Black 30" HR26 1/2 Circle Taper Point 26mm
K833H	D833N	Silk 2-0 Black 30" SH Taperpoint 1/2 Circle 26mm
K834H	D834N	Silk 0 Black 30" SH Taperpoint 1/2 Circle 26mm
U202H	C202N	Chromic Gut 5-0 27" RB-1 Taperpoint 1/2 Circle 17mm
U203H	C203N	Chromic Gut 4-0 27" RB-1 Taperpoint 1/2 Circle 17mm
U204H	C204N	Chromic Gut 3-0 27" RB-1 Taperpoint 1/2 Circle 17mm
VCP944H	M944N	PolySyn FA™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
VR426	M426N	PolySyn FA™ 4-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm
VR490	M490N	Polysyn FA™ 5-0 18" P-1 Precision Reverse Cutting 3/8c 11mm
VR493	M493N	PolySyn FA™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
VR494	M494N	PolySyn FA™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
VR496	M496N	PolySyn FA™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
VR497	M497N	PolySyn FA™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
VR917	M917N	Polysyn FA™ 2-0 Undyed 36" V-34 Taper Cutting 1/2c 36mm
VR935	M442N	Polysyn FA™ 3-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8c 24mm
VR945	M945N	PolySyn FA™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
X424H	E424N	Polyviolene™ 0 Green 30" CT-1 Taperpoint 1/2 Circle 36mm
X425H	E425N	Polyviolene™ 1 Green 30" CT-1 Taperpoint 1/2 Circle 36mm
X519H	E519N	Polyviolene™ 2 Green 30" OS-4 Reverse Cutting 1/2C 22mm
X523H	EX523D	Polyviolene™ 2-0 Green Braided 36" SH Taperpoint 1/2 Circle 26mm
X524H	E524N	Polyviolene™ 0 Green 36" SH Taperpoint 1/2 Circle 26mm Double Armed
X551H	EX551D	Polyviolene™ 4-0 Green Braided 30" RB-1 Taperpoint 1/2 Circle 17mm
X556H	EX556D	Polyviolene™ 5-0 Green Braided 36" RB-1 Taperpoint 1/2 Circle 17mm
X557H	EX557D	Polyviolene™ 4-0 Green Braided 36" RB-1 Taperpoint 1/2 Circle 17mm
X558H	EX558D	Polyviolene™ 3-0 Green Braided 36" RB-1 Taperpoint 1/2 Circle 17mm
X563H	EX563D	Polyviolene™ 2-0 Green Braided 30" SH Taperpoint 1/2 Circle 26mm
X762H	EX762D	Polyviolene™ 3-0 Green Braided 36" SH-1 Taperpoint 1/2 Circle 22mm
X763H	EX763D	Polyviolene™ 2-0 Green Braided 36" SH-1 Taperpoint 1/2 Circle 22mm
X833H	E833N	Polyviolene™ 2-0 Green 30" SH Taperpoint 1/2 Circle 26mm



# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
X834H	E412N	Polyviolene™ 0 Green 30" CT-2 Taperpoint 1/2 Circle 26mm
X834H	E834N	Polyviolene™ 0 Green 30" SH Taperpoint 1/2 Circle 26mm
X927H	EX927D	Polyviolene™ 2-0 Green Braided 30" V-4 Tapercut 3/8 Circle 17mm
X977H	EX977D	Polyviolene™ 2-0 Green Braided 30" V-7 Tapercut 1/2 Circle 26mm
Y213H	Y213D	Monoderm™ 5-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm
Y214H	Y214N	Monoderm™ 4-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm
Y215H	Y215N	Monoderm™ 3-0 Undyed 27" RB-1 Taperpoint 1/2 Circle 17mm
Y266H	Y266N	Monoderm™ 2-0 Undyed 27" CP-1 Reverse Cutting 1/2C 36mm
Y303H	Y303D	Monoderm™ 5-0 Violet 27" RB-1 Taper Point 17mm 1/2C
Y304H	Y304D	Monoderm™ 4-0 Violet 27" RB-1 Taper Point 17mm 1/2C
Y305H	Y305D	Monoderm™ 3-0 Violet 27" RB-1 Taper Point 17mm 1/2C
Y316H	Y316N	Monoderm™ 3-0 Violet 27" SH Taperpoint 1/2 Circle 26mm
Y398H	Y398N	Monoderm™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm
Y399H	Y399N	Monoderm™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm
Y415H	Y415N	Monoderm™ 4-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm
Y416H	Y416N	Monoderm™ 3-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm
Y417H	Y417D	Monoderm™ 2-0 Undyed SH Taper Point 1/2c 26mm
Y417H	Y417N	Monoderm™ 2-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm
Y426H	Y426N	Monoderm™ 4-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm
Y427H	Y427N	Monoderm™ 3-0 Undyed 27" PS-2 Precision Reverse Cutting 3/8C 19mm
Y463G	Y463N	Monoderm™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm
Y464G	Y464N	Monoderm™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm
Y492G	Y492N	Monoderm™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
Y493G	Y493N	Monoderm™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
Y494G	Y494N	Monoderm™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8C 13mm
Y495G	Y495N	Monoderm™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Y496G	Y496N	Monoderm™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Y497G	Y497N	Monoderm™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Y513G	Y513N	Monoderm™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Y682H	Y682N	Monoderm™ 4-0 Undyed 18" PS-1 Precision Reverse Cutting 3/8C 24mm
Y823G	Y823N	Monoderm™ 4-0 Undyed 18" PC-5 Precision Conventional Cutting 3/8C 19mm
Y834G	Y834N	Monoderm™ 5-0 Undyed 18" PC-1 Precision Conventional Cutting 3/8C 13mm
Y844G	Y844N	Monoderm™ 5-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8C 16mm

# Ethicon Cross Reference Guide

Ethicon	SharpPoint Plus™	Product Description
Y845G	Y845N	Monoderm™ 4-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8C 16mm
Y935H	Y935N	Monoderm™ 4-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm
Y936H	Y936N	Monoderm™ 3-0 Undyed 27" PS-1 Precision Reverse Cutting 3/8C 24mm
Y944H	Y944N	Monoderm™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 circle 36mm
Y945H	Y945N	Monoderm™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
Y946H	Y946N	Monoderm™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
Z126H	R126N	PDO™ 5-0 Violet 30" C-1 Taperpoint 3/8 Circle 13mm Double Armed
Z127H	R127N	PDO™ 6-0 Violet 30" C-1 Taperpoint 3/8 Circle 13mm Double Armed

Z148H	R148N	PDO™ 5-0 Violet 30" RB-2 Taperpoint 1/2 Circle 13mm Double Armed
Z259H	R259N	PDO™ 2-0 Undyed 27" CT-1 Taperpoint 1/2 Circle 36mm
Z303H	R303N	PDO™ 5-0 Violet 27" RB-1 Taperpoint 1/2 Circle 17mm
Z304H	R304N	PDO™ 4-0 Violet 27" RB-1 Taperpoint 1/2 Circle 17mm
Z305H	R305N	PDO™ 3-0 Violet 27" RB-1 Taperpoint 1/2 Circle 17mm
Z316H	R316N	PDO™ 3-0 Violet 27" SH Taperpoint 1/2 Circle 26mm
Z320H	R320N	PDO™ 5-0 Violet 30" RB-1 Taperpoint 1/2 Circle 17mm Double Armed
Z333H	R333N	PDO™ 2-0 Violet 27" CT-2 Taperpoint 1/2 Circle 26mm
Z333H	R317N	PDO™ 2-0 Violet 27" SH Taperpoint 1/2 Circle 26mm
Z334H	R334N	PDO™ 0 Violet 27" CT-2 Taperpoint 1/2 Circle 26mm
Z339H	R339N	PDO™ 2-0 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm
Z340H	R340N	PDO™ 0 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm
Z341H	R341N	PDO™ 1 Violet 27" CT-1 Taperpoint 1/2 Circle 36mm
Z346H	R346N	PDO™ 0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
Z347H	R347N	PDO™ 1 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
Z358T	R358N	PDO™ 0 Violet 36" CT Taperpoint 1/2 Circle 40mm
Z359T	R359N	PDO™ 1 Violet 36" CT Taperpoint 1/2 Circle 40mm
Z370T	R370N	PDO™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm
Z371T	R371N	PDO™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm
Z416H	R416N	PDO™ 3-0 Undyed 27" SH Taperpoint 1/2 Circle 26mm
Z422H	R422N	PDO™ 4-0 Violet 27" FS-2 Reverse Cutting 3/8C 19mm
Z423H	R423N	PDO™ 3-0 Violet 27" FS-2 Reverse Cutting 3/8C 19mm
Z443H	R443N	PDO™ 2-0 Violet 27" FS-1 Reverse Cutting 3/8C 24mm
Z493G	R493N	PDO™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm
Z494G	R494N	PDO™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8C 13mm
Z496G	R496N	PDO™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Z497G	R497N	PDO™ 3-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Z513G	R513N	PDO™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8C 19mm
Z990G	R990N	PDO™ 0 Violet 60" CTX Taperpoint 1/2 Circle 48mm

N=box of 12; D=box of 36

# Covidien Cross Reference Guide

Covidien	Sharp Plus™	Product Description
CC923	Y945N	Monoderm™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
CC924	Y946N	Monoderm™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
CC975	Y399N	Monoderm™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm
CL904	G370N	PolySyn™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm
CL905	G371N	PolySyn™ 1 Violet 36" CTX Taperpoint 1/2 Circle 48mm
CL914	G358N	PolySyn™ 0 Violet 36" CT Taperpoint 1/2 Circle 40mm
CL922	G344N	PolySyn™ 3-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
CL923	G345N	PolySyn™ 2-0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
CL924	G346N	PolySyn™ 0 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
CL925	G347N	PolySyn™ 1 Violet 36" CT-1 Taperpoint 1/2 Circle 36mm
CL952	G944N	PolySyn™ 3-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
CL953	G945N	PolySyn™ 2-0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
CL954	G946N	PolySyn™ 0 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
CL955	G947N	PolySyn™ 1 Undyed 36" CT-1 Taperpoint 1/2 Circle 36mm
CL963	G957N	PolySyn™ 2-0 Undyed 36" CT Taperpoint 1/2 Circle 40mm
CL964	G958N	PolySyn™ 0 Undyed 36" CT Taperpoint 1/2 Circle 40mm
CL974	G978N	PolySyn™ 0 Undyed 36" CTX Taperpoint 1/2 Circle 48mm
CL975	G977N	PolySyn™ 1 Undyed 36" CTX Taperpoint 1/2 Circle 48mm
CM904	Y398N	Monoderm™ 0 Violet 36" CTX Taperpoint 1/2 Circle 48mm
CP411	J8411N	Polypropylene 2-0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm
CP412	J8412N	Polypropylene 0 Blue 30" CT-2 Taperpoint 1/2 Circle 26mm
CP412	J8418N	Polypropylene 0 Blue 30" MO-6 Taperpoint 1/2 Circle 26mm
CP423	J8423N	Polypropylene 2-0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm
CP424	J8424N	Polypropylene 0 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm
CP425	J8425N	Polypropylene 1 Blue 30" CT-1 Taperpoint 1/2 Circle 36mm
CP434	J8434N	Polypropylene 0 Blue 30" CT Taperpoint 1/2 Circle 40mm
CP435	J8435N	Polypropylene 1 Blue 30" CT Taperpoint 1/2 Circle 40mm
CP455	J8455N	Polypropylene 1 Blue 30" CTX Taperpoint 1/2 Circle 48mm
CS424	D424N	Silk 0 Black 30" CT-1 Taperpoint 1/2 Circle 36mm
GS65M	D834N	Silk 0 Black 30" SH Taperpoint 1/2 Circle 26mm
GS66M	D833N	Silk 2-0 Black 30" SH Taperpoint 1/2 Circle 26mm
GS67M	D832N	Silk 3-0 Black 30" SH Taperpoint 1/2 Circle 26mm
SC5587G	Y463N	Monoderm™ 5-0 Violet 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SC5590G	Y464N	Monoderm™ 4-0 Violet 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SC5616G	Y495N	Monoderm™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SC5617G	Y496N	Monoderm™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SC5618G	Y497N	Monoderm™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm

# Covidien Cross Reference Guide

Covidien	Sharp point Plus™	Product Description
SC5627G	Y513N	Monoderm™ 4-0 Violet 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SC5688G	Y492N	Monoderm™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SC5689G	Y493N	Monoderm™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SC5690G	Y494N	Monoderm™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SDN5666G	AG667N	Nylon™ 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SL1625G	G489N	PolySyn™ 6-0 Undyed 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm
SL1654	G824N	PolySyn™ 3-0 Undyed 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm
SL1657	G835N	Polysyn™ 4-0 Undyed 18" Precision Conventional Cutting 3/8 16mm
SL1657	G845N	PolySyn™ 4-0 Undyed 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm
SL1659	G843N	Polysyn 5-0 Undyed PC-1 Precision Conventional Cutting 3/8c 13mm
SL5626	G500N	PolySyn™ 5-0 Undyed 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm
SL5627G	G496N	PolySyn™ 4-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SL5628G	G497N	PolySyn™ 3-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SL5632G	G506N	Polysyn 5-0 Undyed 18" PS-4 Precision Reverse Cutting 1/2 Circle 16mm
SL5636G	G495N	PolySyn™ 5-0 Undyed 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SL5687G	G493N	PolySyn™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SL5688G	G492N	PolySyn™ 6-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SL5691G	G494N	PolySyn™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SL5633G	G507N	PolySyn™ 4-0 Undyed 18" PS-4 Precision Reverse Cutting 1/2 Circle 16mm
SL5640G	G683N	PolySyn™ 3-0 Undyed 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm
SN658	A1629N	Nylon™ 4-0 Black 18" FS-1 Reverse Cutting 3/8 Circle 24mm
SN660	A0181N	Nylon™ 6-0 Black 18" FS-3 Reverse Cutting 3/8c 16mm
SN661	AC-0568D	Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8c 19mm
SN661G	A661N	Nylon™ 5-0 Black 18" FS-2 Reverse Cutting 3/8 Circle 19mm
SN662G	A662N	Nylon™ 4-0 Black 18" FS-2 Reverse Cutting 3/8 Circle 19mm
SN663G	A663N	Nylon™ 3-0 Black 18" FS-1 Reverse Cutting 3/8 Circle 24mm
SN669	A669N	Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8 Circle 24mm
SN669	AC-0586D	Nylon™ 3-0 Black 30" FS-1 Reverse Cutting 3/8c 24mm
SN671G	A1671N	Nylon™ 3-0 Black 30" FSL Reverse Cutting 3/8 Circle 30mm
SN671G	A1691N	Nylon™ 3-0 Black 30" PSL Reverse Cutting 3/8 Circle 30mm
SN764G	A664N	Nylon™ 2-0 Black 18" FS Reverse Cutting 3/8 Circle 26mm
SN1693	A1893N	Nylon™ 3-0 Black 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm
SN1964	A1864N	Nylon™ 4-0 Black 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm
SN1994	A1894N	Nylon™ 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm
SN1994	A1994N	Nylon™ (Pliable) 4-0 Black 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm
SN3697	A697N	Nylon™ 6-0 Black 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm
SN3965	A1865N	Nylon™ 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm
SN3965	A1965N	Nylon™ (Pliable) 5-0 Black 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm

N=box of 12; D=box of 36

# Covidien Cross Reference Guide

Covidien	Sharp point Plus™	Product Description
SN5662G	A1662N	Nylon™ 4-0 Black 18" PS-4 Precision Reverse Cutting 1/2 Circle 16mm
SN5663G	A1663N	Nylon™ 3-0 Black 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm
SN5665G	A1281N	Nylon™ 6-0 Black 18" PS-2 Precision Reverse Cutting 3/8c 19mm
SN5666G	A1666N	Nylon™ 5-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SN5666G	AD1666	Nylon™ 5-0 Black 18" PS-2 Precision Reverse Cutting 3/8c 19mm
SN5667	A662D	Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SN5667G	A1667N	Nylon™ 4-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SN5668G	A0556N	Nylon™ 5-0 Black 18" FS-3 Reverse Cutting 3/8c 16mm
SN5668G	A1668N	Nylon™ 5-0 Black 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm
SN5668	A0556N	Nylon™ 5-0 Black 18" FS-3 Reverse Cutting 3/8c 16mm
SN5669G	A1669N	Nylon™ 3-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SN5670	A593N	Nylon™ 2-0 Black 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SN5690G	A690N	Nylon™ 5-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SN5691	A691N	Nylon™ 4-0 Undyed 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SN5696G	A1698N	Nylon™ 6-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SN5698G	A698N	Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SN5698G	AD698N	Nylon™ 5-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
SN5699G	A699N	Nylon™ 4-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SP623	J8623N	Polypropylene 2-0 Blue 30" KS Straight Cutting 60mm
SP623BC	J8623N	Polypropylene 2-0 Blue 30" KS Straight Cutting 60mm
SP661G	J8661N	Polypropylene 5-0 Blue 18" FS-2 Reverse Cutting 3/8 Circle 19mm
SP665G	J8665N	Polypropylene 3-0 Blue 18" FS-2 Reverse Cutting 3/8 Circle 19mm
SP683G	J8683N	Polypropylene 4-0 Blue 18" FS-2 Reverse Cutting 3/8 Circle 19mm
SP684G	J8684N	Polypropylene 3-0 Blue 18" FS-1 Reverse Cutting 3/8 Circle 24mm
SP868G	J8685N	Polypropylene 2-0 Blue 18" FS Reverse Cutting 3/8 Circle 26mm
SP1634	J8634N	Polypropylene 4-0 Blue 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm
SP1635	J8635N	Polypropylene 5-0 Blue 18" PC-3 Precision Conventional Cutting 3/8 Circle 16mm
SP1643	J8632N	Polypropylene 3-0 Blue 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm
SP1644	J8631N	Polypropylene 4-0 Blue 18" PC-5 Precision Conventional Cutting 3/8 Circle 19mm
SP5663G	J8663N	Polypropylene 3-0 Blue 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm
SP5680G	J8680N	Polypropylene 6-0 Blue 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm
SP5681G	J8681N	Polypropylene 5-0 Blue 18" PS-3 Precision Reverse Cutting 3/8 Circle 16mm
SP5682G	J8682N	Polypropylene 4-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SP5686G	J8686N	Polypropylene 5-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SP5687G	J8687N	Polypropylene 3-0 Blue 18" PS-2 Precision Reverse Cutting 3/8 Circle 19mm
SP5695G	J8695N	Polypropylene 6-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm

N=box of 12; D=box of 36

# Covidien Cross Reference Guide

Covidien	Sharpoint Plus™	Product Description
SP5697G	J8697N	Polypropylene 6-0 Blue 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm
SP5698G	J8698N	Polypropylene 5-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SP5699	J8699N	Polypropylene 4-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SP5699G	J8699N	Polypropylene 4-0 Blue 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SS623	D623D	Silk 2-0 Black 1x30" KS Straight Cutting 60mm
SS678	D580N	Silk 0 Black 18" PSL Reverse Cutting 3/8 Circle 30mm
SS678	D678N	Silk 0 Black 18" FSL Reverse Cutting 3/8 Circle 30mm
SS679	D679N	Silk 2-0 Black 30" FSL Reverse Cutting 3/8c 30mm
SS684G	D684N	Silk 3-0 Black 18" FS-1 Reverse Cutting 3/8 Circle 24mm
SS685G	D685N	Silk 2-0 Black 18" FS Reverse Cutting 3/8 Circle 26mm
SS70M	D590N	Silk 0 Black 30" PSL Reverse Cutting 3/8 Circle 30mm
SS70M	D680N	Silk 0 Black 30" FSL Reverse Cutting 3/8 Circle 30mm
SS1639G	D639N	Silk 6-0 Black 18" P-1 Precision Reverse Cutting 3/8 Circle 11mm
SS1639G	D1639N	Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
SS5639	D1639N	Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
SS5640	D640N	Silk 5-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SS5641G	D641N	Silk 4-0 Black 18" P-3 Precision Reverse Cutting 3/8 Circle 13mm
SS5649G	D1639N	Silk 6-0 Black 18" P-3 Precision Reverse Cutting 3/8c 13mm
SS5677	D683N	Silk 4-0 Black 1x18" PS-2 Precision Reverse Cutting 3/8c 19mm
SS5679	D1679N	Silk 3-0 Black 18" FS-2 Reverse Cutting 3/8c 19mm
SS5684	D1684N	Silk 3-0 Black 18" PS-1 Precision Reverse Cutting 3/8 Circle 24mm
SS5685G	D0238N	Silk 2-0 Black 18" FS-1 Reverse Cutting 3/8c 24mm
SV9915	M490N	Polysyn™ FA 5-0 18" P-1 Precision Reverse Cutting 3/8c 11mm
VP37MX	J8755N	Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP72MX	J8747N	Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP74X	J8610N	Polypropylene 6-0 Blue 24" BV Taperpoint 3/8 Circle 11mm Double Armed
VP205X	J8554D	Polypropylene 5-0 Blue 24" RB-2 Taperpoint 1/2 Circle 13mm Double Armed
VP422X	J8522D	Polypropylene 3-0 Blue 36" SH Taperpoint 1/2 Circle 26mm Double Armed
VP521X	J8521D	Polypropylene 4-0 Blue 36" SH Taperpoint 1/2 Circle 26mm Double Armed
VP523	J8523D	Polypropylene 2-0 Blue 36" SH Taperpoint 1/2 Circle 26mm Double Armed
VP542X	J8534D	Polypropylene 3-0 Blue 48" SH Taperpoint 1/2 Circle 26mm Double Armed
VP543	J8533D	Polypropylene 2-0 Blue 48" SH Taperpoint 1/2 Circle 26mm Double Armed
VPF543	J8533D	Polypropylene 2-0 Blue 48" SH Taperpoint 1/2 Circle 26mm Double Armed
VP556X	J8556D	Polypropylene 5-0 Blue 36" RB-1 Taperpoint 1/2 Circle 17mm Double Armed

N=box of 12; D=box of 36

# Covidien Cross Reference Guide

Covidien	Sharp Plus™	Product Description
VP557X	J8557D	Polypropylene 4-0 Blue 36" RB-1 Taperpoint 1/2 Circle 17mm Double Armed
VP558X	J8558D	Polypropylene 3-0 Blue 36" RB-1 Taperpoint 1/2 Circle 17mm Double Armed
VP566X	J8555D	Polypropylene 5-0 24" RB-1 Taper Point 1/2 Circle 17mm
VP580X	J8580D	Polypropylene 5-0 Blue 36" BB Taperpoint 3/8 Circle 17mm Double Armed
VP581X	J8581D	Polypropylene 4-0 Blue 36" BB Taperpoint 3/8 Circle 17mm Double Armed
VP701X	J8701D	Polypropylene 7-0 Blue 18" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP702X	J8702D	Polypropylene 7-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP703X	J8805D	Polypropylene 6-0 Blue 24" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP706X	J8706D	Polypropylene 6-0 Blue 30" C-1 Taperpoint 3/8 Circle 13mm Double Armed
VP709X	J8709D	Polypropylene 6-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP710X	J8710D	Polypropylene 5-0 Blue 30" RB-2 Taperpoint 1/2 Circle 13mm Double Armed
VP713X	J8703D	Polypropylene 7-0 Blue 30" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP714X	J8714D	Polypropylene 6-0 18" RB-2 Taper Point 1/2 Circle 13mm
VP715X	J8713D	Polypropylene 5-0 18" RB-2 Taper Point 1/2 Circle 13mm
VP718X	J8718D	Polypropylene 6-0 Blue 18" C-1 Taperpoint 3/8 Circle 13mm Double Armed
VP719X	J8717D	Polypropylene 5-0 Blue 18" C-1 Taperpoint 3/8 Circle 13mm Double Armed
VP720X	J8720D	Polypropylene 5-0 Blue 36" C-1 Taperpoint 3/8 Circle 13mm Double Armed
VP725X	J8725D	Polypropylene 5-0 Blue 24" C-1 Taperpoint 3/8 Circle 13mm Double Armed
VP726X	J8726D	Polypropylene 6-0 Blue 24" C-1 Taperpoint 3/8 Circle 13mm Double Armed
VP733X	J8711D	Polypropylene 6-0 Blue 30" RB-2 Taperpoint 1/2 Circle 13mm Double Armed
VP746X	J8806D	Polypropylene 6-0 Blue 18" BV-1 Taperpoint 3/8 Circle 9.3mm Double Armed
VP760X	J8762D	Polypropylene 3-0 Blue 30" SH-1 Taperpoint 1/2 Circle 22mm Double Armed
VP776X	J8776N	Polypropylene 6-0 Blue 30" BV Taperpoint 3/8 Circle 11mm Double Armed
VP832X	J8832N	Polypropylene 3-0 Blue 30" SH Taperpoint 1/2 Circle 26mm
VP833	J8833N	Polypropylene 2-0 Blue 30" SH Taperpoint 1/2 Circle 26mm
VP834	J8418N	Polypropylene 0 Blue 30" MO-6 Taperpoint 1/2 Circle 26mm
VP834	J8834N	Polypropylene 0 Blue 30" SH Taperpoint 1/2 Circle 26mm
VP842X	J8842D	Polypropylene 3-0 Blue 36" MH Taperpoint 1/2 Circle 36mm Double Armed
VP843	J8843D	Polypropylene 2-0 Blue 36" MH Taperpoint 1/2 Circle 36mm Double Armed
VP853	J8853D	Polypropylene 2-0 48" MH Taper Point 1/2 Circle 36mm













Order Online [www.SurgicalSpecialties.com](http://www.SurgicalSpecialties.com)



Order by Phone 877-991-1110

---

SharpPoint™ is also available through major healthcare distributors.

**SURGICAL  
SPECIALTIES**  
CORPORATION™

For More Information, Visit [SurgicalSpecialties.com](http://SurgicalSpecialties.com)

Surgical Specialties Corporation, 247 Station Drive, Suite NE1, Westwood, MA 02090, U.S.A.  
Phone: 610.404.3520 • 877.991.1110 (USA and Canada only) • Fax: 610.404.4010  
All trademarks are property of their respective owners. ©2020 Surgical Specialties Corporation.  
All Rights Reserved. SSC-028 R8 01/2021